



TruPass™ SW (Short-Wire) Triple Lumen Papillotome

Compatible with short-wire ERCP techniques, CONMED's NEW TruPass™ SW Papillotomes are designed to deliver guaranteed orientation and currently have one of the smallest tips of any .035" guidewire-compatible triple lumen papillotomes available.

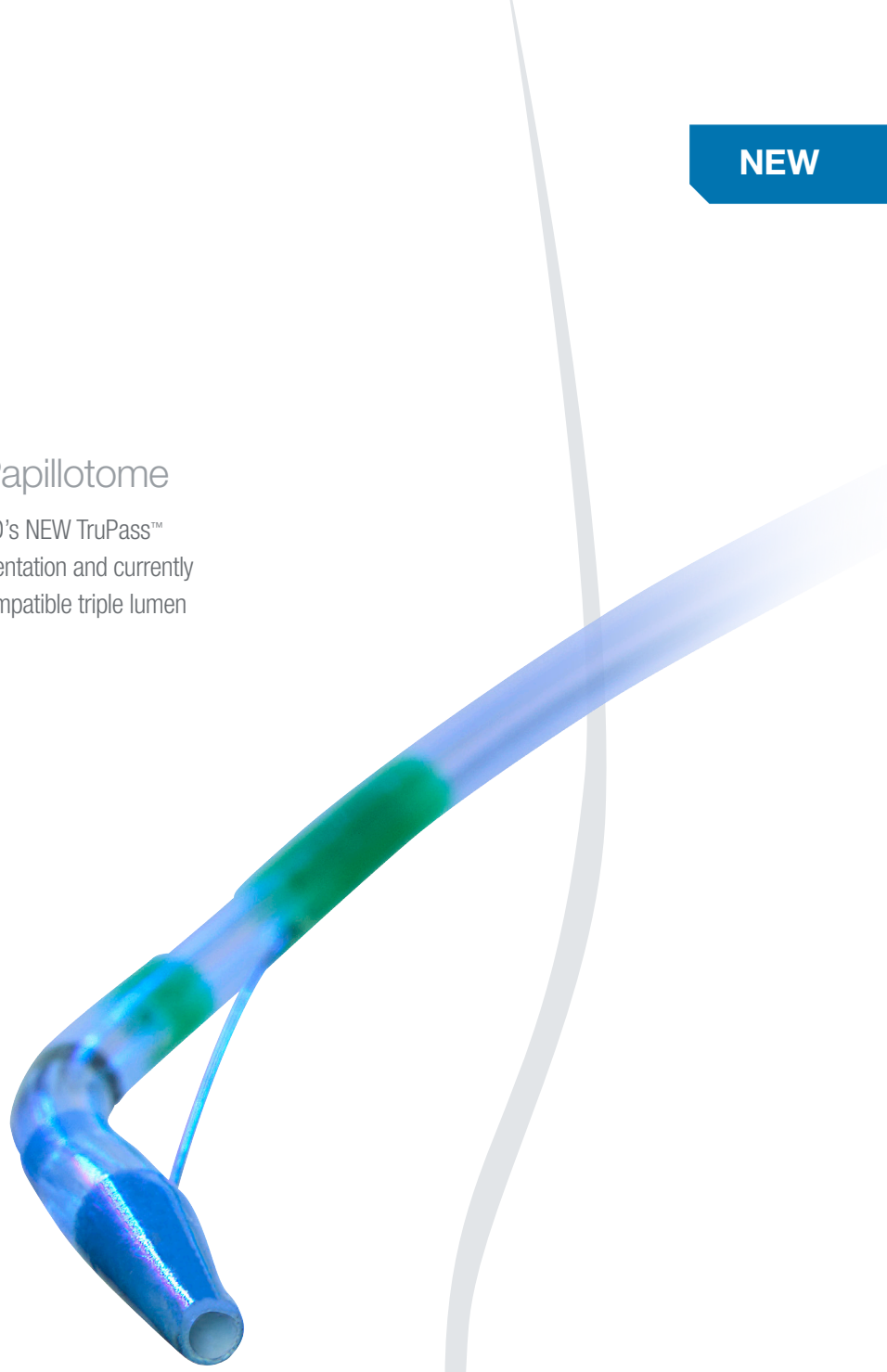


GORE® VIABIL®
BILIARY ENDOPROSTHESIS

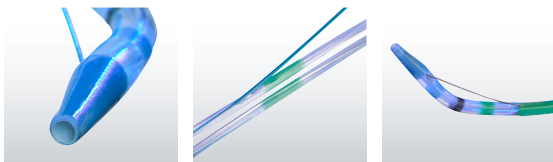

ACE™
BILIARY GUIDEWIRE



FLEXXUS®
ENDOSCOPIC BILIARY STENT

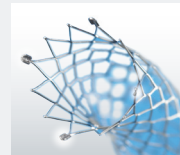

To learn more about these and other innovative products, call **800-448-6506** or visit **CONMED.com**.

Guaranteed Orientation: TruPass™ Papillotomes are guaranteed to orient between 11 and 1 o'clock to help simplify ERCP procedures.

Smaller Tip: Currently one of the smallest .035" guidewire-compatible triple lumen papillotomes available, the 4.3Fr tip is designed to improve cannulation performance and maneuverability.

Company	Product	Distal Tip Diameter (Fr)
CONMED	TruPass™	4.3
Boston Scientific	Autotome™	4.4*
Olympus	CleverCut3V	4.4**

Integrated Lumen Technology for Optimal Guidewire Positioning: To facilitate cannulation, the guidewire exits the distal tip of the TruPass™ sphinctertome through a central lumen. This helps to position the guidewire axis more preferentially to help with guidewire cannulation.

Color-Coded Cutting Wire Markers: Easy visibility for determining papillotome depth and cutting wire location.

Smooth, Atraumatic Radiopaque Tip: The smooth tapered tip reduces tissue trauma to the papilla. Tips are also radiopaque to enable fluoroscopic visualization.

▶ **For a complete ERCP Solution, use TruPass™ SW Triple Lumen Papillotomes with FIT™ SW Stone Balloons and FX Wire® Advanced Measurement Guidewire.**

Ordering Information

Description	Cath Length	Tip O.D.	Cath O.D.	Cutting Wire Length	Cutting Wire Type	Guidewire Compatibility	Qty	Catalog Number
NEW TruPass™ Triple Lumen Papillotome								
TruPass™ Short Nose, 8mm	190cm	4.3 F	7-4.3 F	20mm	Monofilament	035" / .89mm	1/box	T2043F
TruPass™ Short Nose, 8mm	190cm	4.3 F	7-4.3 F	30mm	Monofilament	035" / .89mm	1/box	T3043F

Color Coded Guidewire Port, Exit Marker Band, Cutting Wire Markers.

To learn more about this and other innovative products, call **1-800-448-6506** or visit CONMED.com.

Customer Service: 1-800-448-6506
International Sales: 1-315-797-8375
Fax: 1-800-438-3051

CONMED.com
info@mail.conmed.com

* Boston Scientific Autome™ RX brochure. ENDO-377606-AA May 2016
** Olympus, Medical.olympusamerica.com/products/clevercut-3v-digital-wire-guided, Accessed on April 20, 2017