




ANCHOR^{*}

Tissue Retrieval System™

New Sizes
Offer **40%** More
Volume!

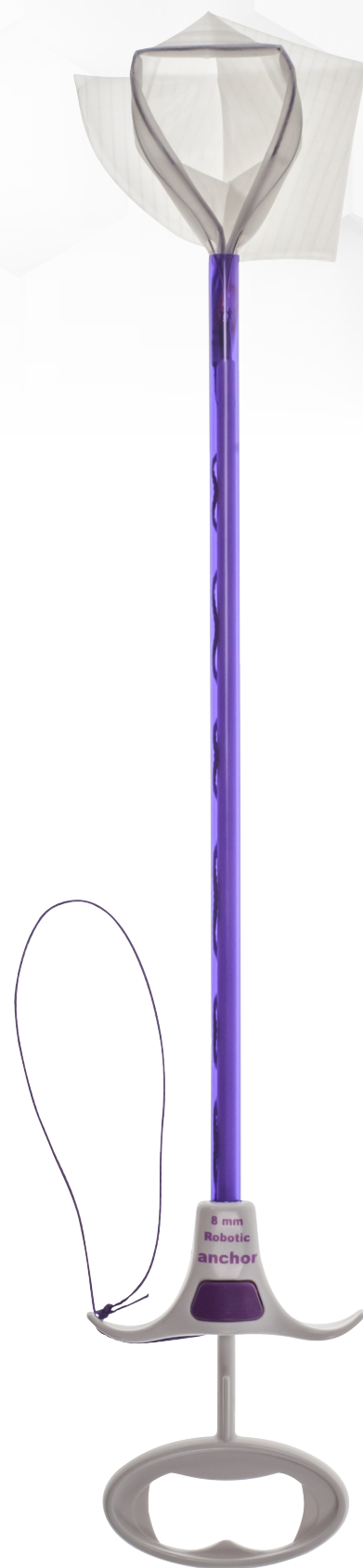

Ready to ask for the “purple one”?

Contact your local CONMED Representative
for a demonstration or visit www.CONMED.com for more information.



6 BAG VOLUME OPTIONS

125 ml • 175 ml • 250 ml • 1200 ml • 1800 ml • 6000 ml

TRS-ROBO-8
125 mL



TRS-ROBO-8L
175 mL



TRS100SB2
175 mL



TRS100SB2L
250 mL



TRS-ROBO-12
250 mL



TRS-VATS-15
1200 mL



TRS175SB2
1200 mL



TRS190SB2
1800 mL



TRS-TV-25
6000 mL



FEATURES & BENEFITS

New Sizes!

- New 8mm and 10mm options feature 40% more bag volume, accommodating a wider range of specimens.

Rip-Stop Nylon

- Proprietary material provides exceptional strength and quality, preventing breaks, leaks, and tears.
- Bag shapes specifically designed to accommodate a wide range of procedures and specialties (thoracic, gynecology, urology, etc.)

Unique Deployment Mechanism

- Upon deployment, Anchor automatically opens and unrolls, offering increased procedural efficiency.
- Can be redeployed in the same case to retrieve multiple specimens, reducing the need for multiple bags.

Simple, Easy-to-Use Design

- Intuitive design provides ease of use during initial, repeat, and final deployment.
- Insert indicator ensures bag will be in proper position when deployed.
- For final deployment, simply **Push-Push-Pull**:
 - **Push** and hold the purple button
 - **Push** the plunger into the introducer until a 'click' is heard
 - **Pull** back on the plunger to remove it from the introducer tube, allowing for final specimen removal.

ORDERING & PRODUCT INFORMATION

TROCAR SIZE	MOUTH DIAMETER	BAG LENGTH	VOLUME	QTY	CATALOG NUMBER
8mm	2in (5.1cm)	5.5in (14cm)	125mL	5/Box	TRS-ROBO-8
8mm (NEW)	2.2in (5.6cm)	5.5in (14cm)	175mL	5/Box	TRS-ROBO-8L
10mm	2.2in (5.6cm)	5.5in (14cm)	175mL	5/Box	TRS100SB2
10mm (NEW)	3in (7.6cm)	5.5in (14cm)	250mL	5/Box	TRS100SB2L
12mm	3in (7.6cm)	5.5in (14cm)	250mL	5/Box	TRS-ROBO-12
15mm	5.5in (14cm)	9.9in (25.4cm)	1200mL	3/Box	TRS-VATS-15
15mm	5.5in (14cm)	9.9in (25.4cm)	1200mL	3/Box	TRS175SB2
15mm	6in (15.2cm)	10.4in (26.5cm)	1800mL	3/Box	TRS190SB2
25mm	6.5in (16.5cm)	15in (38.1cm)	6000mL	3/Box	TRS-TV-25
Specimen Bag Alone					
10mm	2.2in (5.6cm)	5.5in (14cm)	175mL	10/Box	TRS070
10mm	3.5in (8.9cm)	5.5in (14cm)	250mL	10/Box	TRS080
15mm	5.5in (14cm)	9.9in (25.4cm)	1200mL	5/Box	TRS090

* Anchor is a trademark of the Anchor Products Company, Addison, Illinois.

CONMED Corporation
11311 Concept Blvd.
Largo, FL 33773

Toll Free: 1-866-4CONMED
International: 727-214-3000
www.CONMED.com

customerexperience@conmed.com
internationalorders@conmed.com