



Looking Glass[®]


CAMERA HEAD

Customizable, autoclavable, and adjustable from the sterile field — CONMED's Looking Glass[®] camera heads are designed with innovative technology for maximum performance and efficiency.



Features and Benefits

Crisp 4K Technology	4K delivers precise, detailed imaging for multi-specialty visualization. Two unique configurations for your surgeon's preference.
Innovative 3-CMOS Color Imaging	3-CMOS offers a dedicated chip for each color spectrum (red, green, blue), while a single CMOS chip works to blend all colors.
Customizable 5-Button Design	Unique 5-button design allows customization of the system for easy preferences and adjustments from the OR sterile field.
Unique ShockFlex[™] Durability	Resilient ShockFlex [™] prism mounting prevents damage from accidental drops.
Autoclavable to Keep Costs Down	Autoclavable camera heads eliminate the need for expensive sterilization methods. This improves efficiency while maintaining maximum performance.
7-Year Warranty Against Autoclave Failure	Provides peace of mind against steam penetration to ensure durability of your purchase

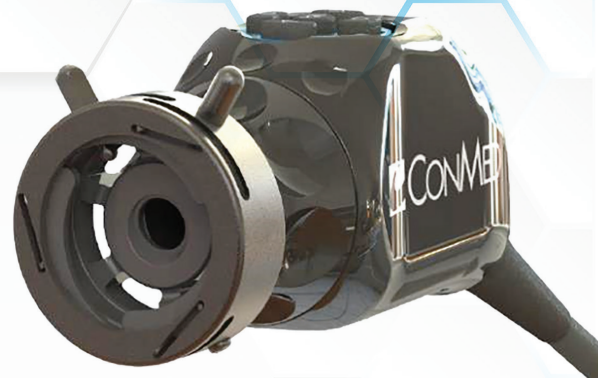

Product Specifications

CMOS Sensors	3 CMOS sensors 1/3 inch
Electronic Zoom	Up to 2:1, controllable by camera head button
IPX Rating	IPX7
Sterilization	See IFU for sterilization instructions ²

Ordering Information

Description	Catalog Number
Looking Glass® Camera Head — 4K 3MOS Eyecup	AV41CHE
Looking Glass® Camera Head — 4K 3MOS Eyecup — Locking	AV41CHE-L

Learn how you can invest in best-in-class Looking Glass® innovation and efficiency by scheduling an evaluation with your local CONMED Sales Representative today at **1-866-4CONMED** or **www.CONMED.com**.

CONMED Corporation
11311 Concept Blvd.
Largo, FL 33773

Toll Free: 1-866-4CONMED
International: 727-214-3000
www.CONMED.com

customerexperience@conmed.com
internationalorders@conmed.com