




NEW

## Lanza™ Tapered Specialty Shaver Blade

CONMED's NEW Lanza™ Tapered Specialty Shaver Blade has a tapered distal tip, smooth cutting profile and dual-inner cutting window that is designed to provide precise, efficient resection in tight anatomical spaces and sensitive areas.

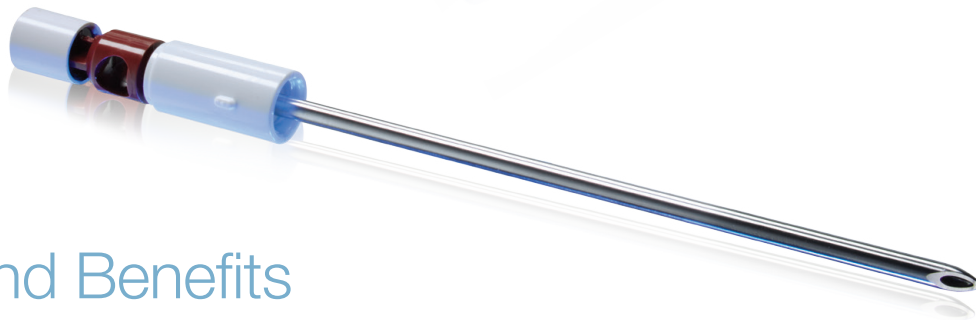

## Features and Benefits

### Better Access

**Tapered Distal Tip** – Designed to provide easy access to tight spaces, difficult-to-reach anatomy and around sensitive areas. Lanza can also serve as an effective soft tissue probe.

### Faster Cutting

**Dual Inner Cutting Window** – Compared to a standard blade<sup>1</sup>, the dual inner cutting window provides twice as many cuts per rotation to help maximize volumetric resection.

### Reduced Particulate


**TrimRite® Steel & Shaft Lubrication** – Designed to reduce metal particulate generation, the shrink tubing on the inner blade shaft provides a barrier to prevent metal-on-metal contact<sup>2</sup> while the exterior features hardened TrimRite® Stainless Steel that helps maintain a sharp cutting edge and blade integrity.<sup>3</sup>

### Smooth Cutting

**Smooth Cutting Profile** – Engineered to prevent snagging on soft tissue when navigating tight spaces, while the shallow exterior window and dual inner cutting window are designed to reduce clogging.

To learn more about this and other innovative products,  
call **1-800-237-0169** or visit **CONMED.com**.

<sup>1</sup> Standard single cutting window blade <sup>2</sup> TR13-200 <sup>3</sup> TR17-00566-1

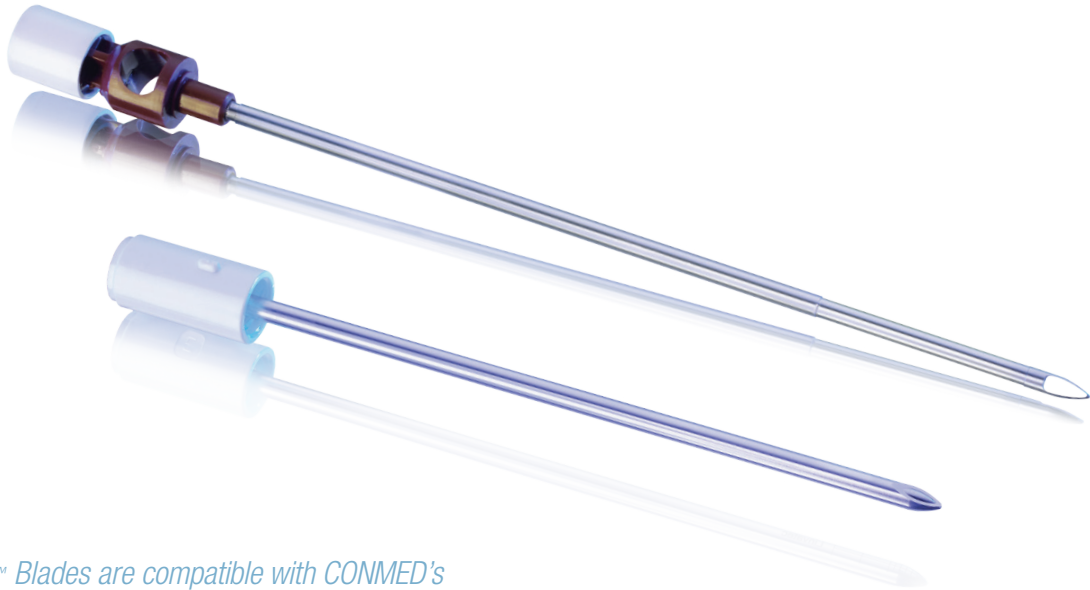

*Lanza™ Blades are compatible with CONMED's Advantage® Turbo and Ergo™ Shaver Handpieces and our D3000 and D4000 consoles.*

## Ordering Information

### Description

### Catalog Number

NEW

Lanza® Tapered Specialty Shaver Blade

9391A

#### Advantage® Turbo Shaver Handpiece

Two-Button Advantage® Turbo Shaver

D9924

No-Button Advantage® Turbo Shaver (Footswitch Only)

D9920

#### Ergo™ Shaver Handpiece

Ergo™ Shaver Handpiece, No-Button

D4200

Ergo™ Shaver Handpiece, Two-Button

D4240

#### Electric Consoles and Accessories

D3000 Arthro Power Console

D3000

D4000 Arthro Power Console

D4000

CONMED Corporation  
525 French Road, Utica, NY 13502

Local 727-392-6464  
Toll Free 1-800-237-0169

CONMED.com  
customer\_service@CONMED.com