

OPERATING BEYOND THE EXPECTED™



THE DETACHATIP® II SYSTEM

THE MULTI-USE ENDOSURGERY
INSTRUMENT SOLUTION

 CONMED™
ENDOSURGERY

An Unmatched Combination of Efficiencies

The DetachaTip®II System from ConMed Corporation provides the ultimate balance between a reusable and disposable hand instrument system. This unique two piece handle and integrated shaft instrument concept delivers excellent cleanability, cost reduction, and consistent performance by enabling you to easily replace high quality, multi-use components as needed. Designed to withstand the rigors of multiple use, these instruments dramatically reduce procedural costs. DetachaTip®II eliminates costly repairs and inconsistent performance of conventional reusables. Only DetachaTip®II offers you this combination of efficiencies in performance and affordability.



Ratchet Handle

Compact Ratchet Handle

In-Line Ratchet Handle



Luer Lock Cap

Rotator

accommodates entire fingertip for comfortable, precise control of tip rotation

New Luer Port

permits high-flow cleaning of the central shaft

Advanced Locking Mechanism

prevents unintentional rotation of tip



Meeker Dissector



Tapered Dissector



Curved Dissector



5mm Babcock Grasper




10mm Babcock Grasper



Endow Grasper

Performs to the Highest Standards

- Simple two piece, high quality, multi-use components easily assemble and disassemble with no extra parts
- Convenient, Individual components packaged sterile for immediate introduction into the operative field
- Integrated shaft/tip assembly delivers consistent, smooth performance
- Durable Luer cleaning port permits high-flow flushing of the central shaft for excellent cleaning results
- Comfortable, durable, ergonomically designed handles provides long-term, repair-free service
- Eliminates cost, inconvenience and unreliability of reusable instrument repairs

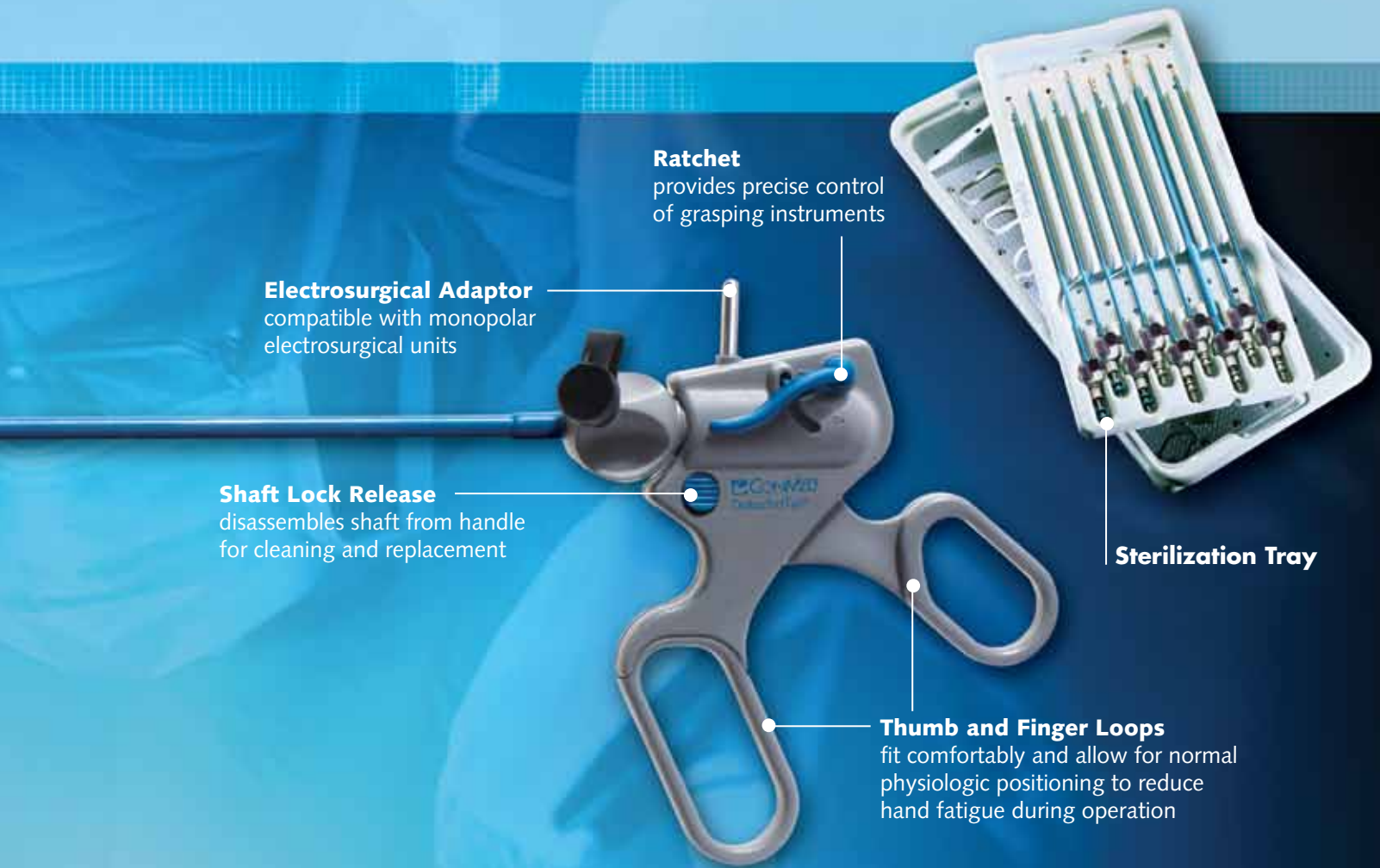

Ratchet
provides precise control
of grasping instruments

Electrosurgical Adaptor
compatible with monopolar
electrosurgical units

Shaft Lock Release
disassembles shaft from handle
for cleaning and replacement

Sterilization Tray

Thumb and Finger Loops
fit comfortably and allow for normal
physiologic positioning to reduce
hand fatigue during operation



leave™
sper

**Allis
Grasper**

**Straight
Fenestrated
Grasper**

**Hook
Scissor**

**Mini-
Metzenbaum
Scissor**

**Metzenbaum
Scissor**

DetachaTip® II

ORDERING INFORMATION



	REFERENCE NO.	PRODUCT DESCRIPTION	PACKAGING
DetachaTip® II Scissors			
	2-1003	Curved Metzenbaum, 33cm x 5mm	1/cs
	2-4301	Curved Metzenbaum, 43cm x 5mm	1/cs
	2-1004	Curved Mini-Metzenbaum, 33cm x 5mm	1/cs
	2-4304	Curved Mini-Metzenbaum, 43cm x 5mm	1/cs
	2-1013	Hook, 33cm x 5mm	1/cs
DetachaTip® II Graspers			
	2-1008	Straight Fenestrated Grasper, 33cm x 5mm	1/cs
	2-4307	Straight Fenestrated Grasper, 43cm x 5mm	1/cs
	2-1028	Endoweave™ Grasper, 33cm x 5mm	1/cs
	2-4328	Endoweave™ Grasper, 43cm x 5mm	1/cs
	2-1019	Allis Grasper, 33cm x 5mm	1/cs
	2-4319	Allis Grasper, 43cm x 5mm	1/cs
	2-1005	5mm Babcock Grasper, 33cm x 5mm	1/cs
	2-4305	5mm Babcock Grasper, 43cm x 5mm	1/cs
	2-1014	10mm Babcock Grasper, 33cm x 10mm	1/cs
DetachaTip® II Dissectors			
	2-1009	Maryland (Curved) Dissector, 33cm x 5mm	1/cs
	2-4308	Maryland (Curved) Dissector, 43cm x 5mm	1/cs
	2-1017	Meeker (Right Angle) Dissector, 33cm x 5mm	1/cs
	2-4317	Meeker (Right Angle) Dissector, 43cm x 5mm	1/cs
	2-1018	Tapered (Dolphin) Dissector, 33cm x 5mm	1/cs
	2-4318	Tapered (Dolphin) Dissector, 43cm x 5mm	1/cs
DetachaTip® Handles			
	1-1010	Standard Handle with Ratchet	1/cs
	1-1015	Compact Handle with Ratchet	1/cs
	1-1024	In-Line Handle with Ratchet	1/cs
DetachaTip® II Accessories			
	1-1027	Three-Part Sterilization Tray, for 33cm DetachaTips	1/cs
	1-4327	Three-Part Sterilization Tray, for 43cm DetachaTips	1/cs
DetachaTip® II			
	CD203	Luer lock Cap	5/Bx

Other reusable products from ConMed Endosurgery™: OnePort® Trocar System and Rotating Head Skin Stapler

CONMED CORPORATION PRODUCT AREAS:

ARTHROSCOPY • ELECTROSURGERY • ENDOSCOPY • ENDOSURGERY • GASTROENTEROLOGY • INTEGRATED SYSTEMS • PATIENT CARE • POWERED INSTRUMENTS • PULMONOLOGY



525 French Road
 Utica, NY 13502
 (315) 797-8375
 Customer Service: 1-800-448-6506
 FAX: 1-800-438-3051
 International Sales: +1(315) 797-8375 • FAX: +1(315) 735-6235
 email: info@mail.conmed.com
 www.conmed.com