


ConMed Endoscope Cleaning Brushes

Channel Master™ and Blue Bullet® Disposable Cleaning Brushes

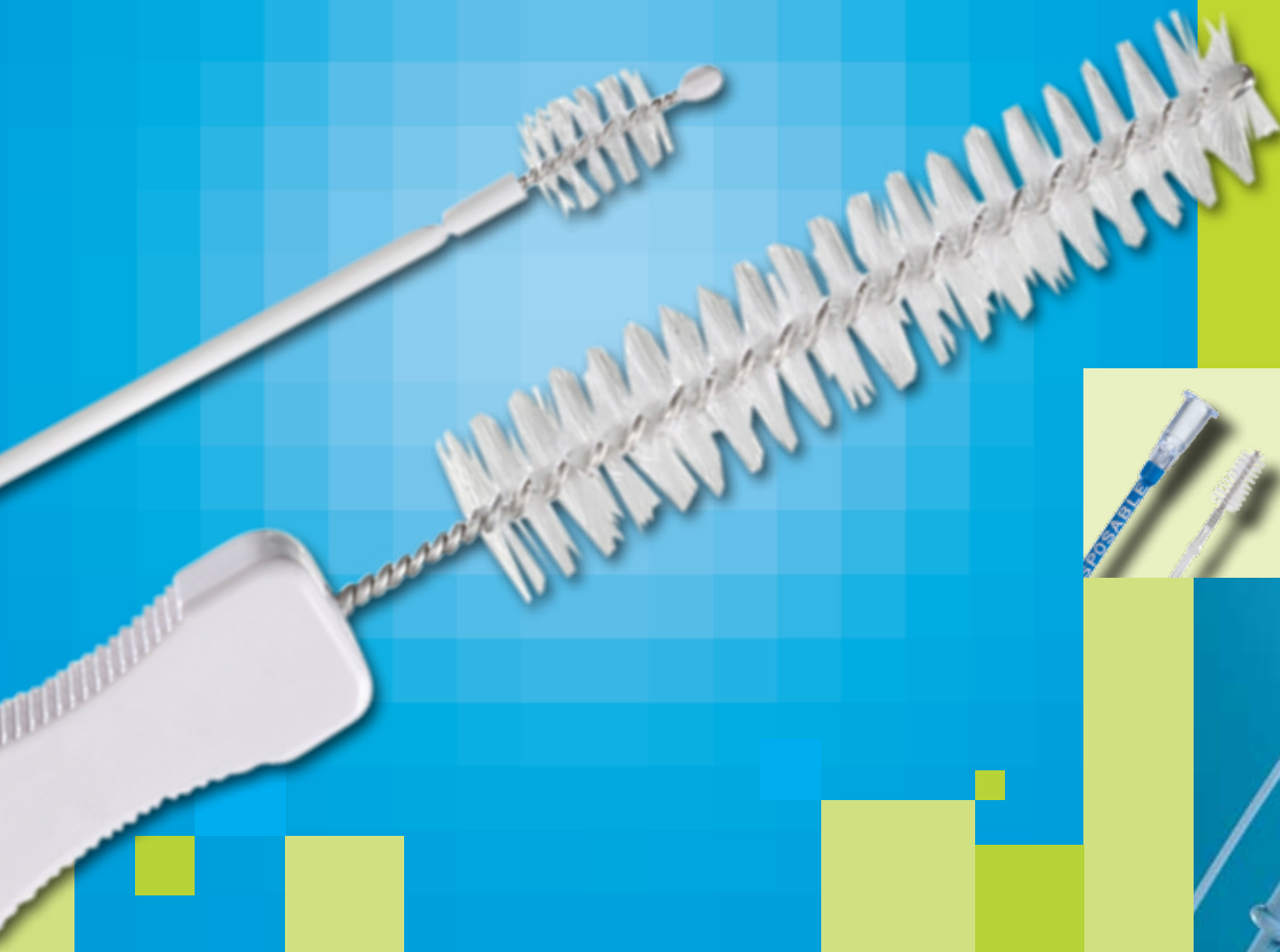


Universal Brush Designed for Optimal Cleaning

Applied Ingenuity™

CONMED™
ENDOSCOPIC TECHNOLOGIES

ConMed Endoscope Cleaning Brushes

- **Universal Brush Design**
 - Cleans full range of endoscope channels
- **Channel Master Custom Crimped Shaft**
 - Effectively traverse tight channel bends
 - Deflects through 90° valve port transitions
- **Rounded Brush Tip**
 - Enhances passage through scope channel
 - Protects biopsy channel from damage
- **Combination Brush**
 - Tapered bristles to clean valve ports, seats and buttons
 - Convenient "no slip" ridges
- **Dual-ended Brush**
 - Utilize two brushes with one catheter pass
- **Blue Bullet® Series**
 - Longer working shaft, and stiffer bristles to effectively clean all endoscope channels

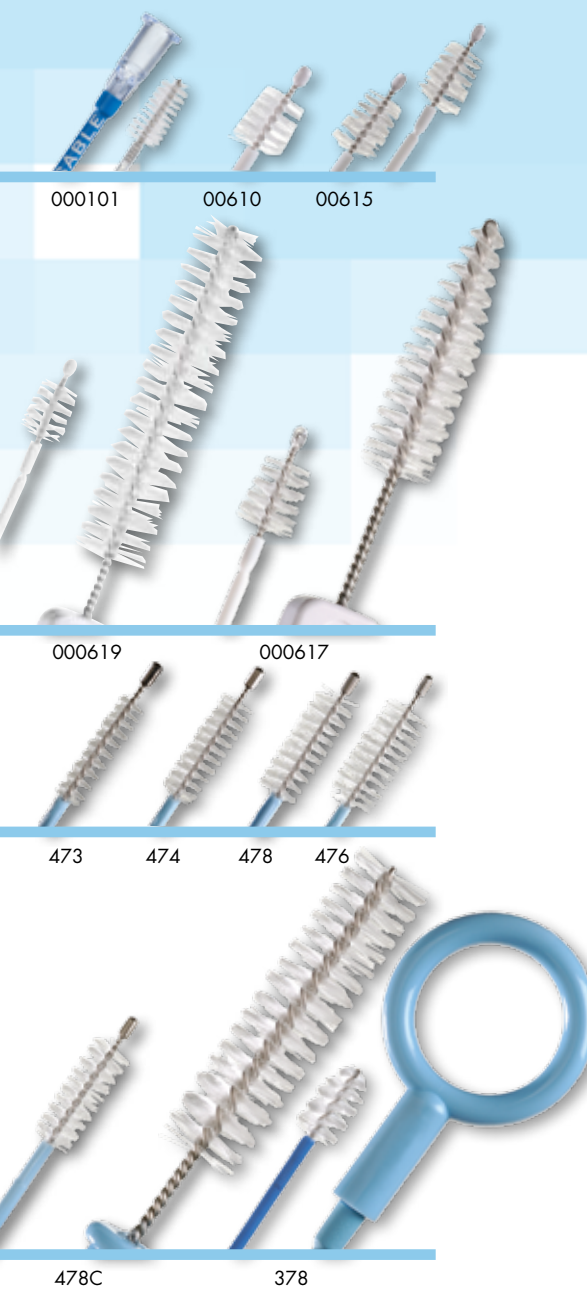

000101 00610 00615

000619 000617

473 474 478 476

478C 378

Channel Master™ Endoscope Cleaning Brushes

Latex **FREE**

Non-Sterile, Disposable - Biopsy Channel Range 2.0-2.4mm, Shaft O.D. 1.7mm

Catalog Number	Scope	Diameter Brush	Length	Length	Shaft	Unit
000610	Pediatric Brush	6.4mm	7.4mm	160cm	50/box	
000615	Dual-Ended Brush	6.4mm	7.4mm	230cm	50/box	
000639	Single-Ended Brush	6.4mm	7.4mm	230cm	50/box	
000617	Combination Valve/Cleaning Brush for Gastroscope/Colonoscope	7.9mm/ 6.4mm	29.2mm/ 7.4mm	230cm	50/box	
000619*	Large Combination Valve/Cleaning Brush for Gastroscope/Colonoscope	9.7mm/ 6.4mm	45.2mm/ 7.4 mm	230cm	50/box	

* Large combination brush includes longer & larger diameter bristles on valve brush.

Whistle® Endoscope Cleaning Brushes

Latex **FREE**

Non-Sterile, Disposable - Biopsy Channel Minimum 2.8mm, Shaft O.D. 1.8mm

Catalog Number	Scope	Diameter Brush	Length	Length	Shaft	Unit
000101	Whistle® Cleaning Brush	6.4mm	13.1mm	230cm	50/box	

Contains female luer adaptor, hollow shaft, and port hole for fluid injection to soften channel obstructions.

Blue Bullet® Disposable Cleaning Brushes

Latex **FREE**

Non-Sterile, Disposable

Catalog Number	Scope	Diameter Brush	Length	Length	Shaft	Unit
473	Gastro/Duodeno/Colonoscope	3mm	13mm	250cm	50/box	
474	Gastro/Duodeno/Colonoscope	4mm	13mm	250cm	50/box	
478	Gastro/Duodeno/Colonoscope	5mm	13mm	250cm	50/box	
476	Gastro/Duodeno/Colonoscope	6mm	13mm	250cm	50/box	
478C	Combination Valve/Channel Cleaning Brush for Gastroscope/Colonoscope	10mm	38mm	250cm	50/box	
478D	Double End Cleaning Brush for Gastroscope/Colonoscope	5mm	13mm	250cm	50/box	

The 400 series Blue Bullet® feature a blue opaque tubing.

Reusable Cleaning Brush

Latex **FREE**

Non-Sterile, Reusable

Catalog Number	Scope	Diameter Brush	Length	Length	Shaft	Unit
378	Gastro/Duodeno/Colonoscope	6mm	9.5mm	230cm	1/box	

The reusable scope cleaning brush features a controllable ring handle and a Teflon® coated spring sheath.

CONMED CORPORATION PRODUCT AREAS:

ARTHROSCOPY • ELECTROSURGERY • ENDOSCOPY • ENDOSURGERY • GASTROENTEROLOGY • INTEGRATED SYSTEMS • PATIENT CARE • POWERED INSTRUMENTS • PULMONOLOGY



Subsidiary of ConMed Corporation

525 French Road
Utica, NY 13502
(315) 797-8375

Customer Service: 1-800-448-6506

FAX: 1-800-438-3051

International Sales: +1(315) 797-8375 • FAX: +1(315) 735-6235

email: info@mail.conmed.com

www.conmed.com

© 2009 ConMed Corporation, 4/2009, Control #12616 Rev A

MBR 9132