


**BEAMER***
ARGON SNARE PROBE

Snare and
Resect It.



Beam It.

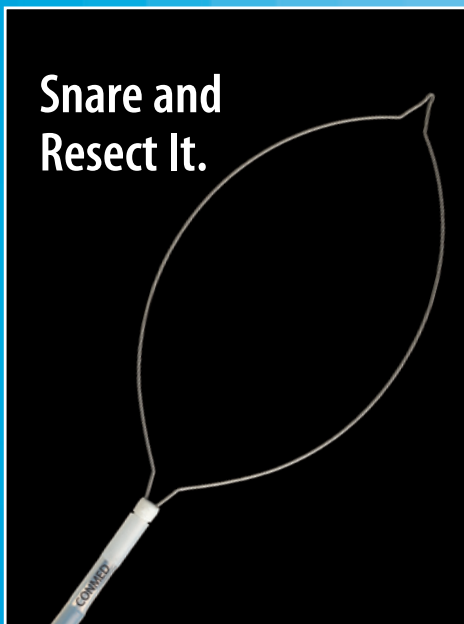

All-in-One.

First in
Argon Beam
Coagulation



A Winning Combination

Applied Ingenuity™

 **CONMED**
ENDOSCOPIC TECHNOLOGIES

BEAMER^{*} ARGON SNARE PROBE

COMBINATION THERAPY

- Interchangeable snare and argon probe allows en bloc or piecemeal resection of lesions followed by argon coagulation for clean up of residual tissue

VERSATILE THERAPY

- Probe beam length, width, and effect can be adjusted by changing wattage, argon flow rate and current mode on Beamer System


CERAMIC TIP

- Guides argon beam flow for targeted therapy
- Offers visual contrast of probe against mucosa

LONG BEAM

- Argon beam can be initiated without contacting the mucosa

BEAMER ARGON SNARE PROBE KEY FEATURES

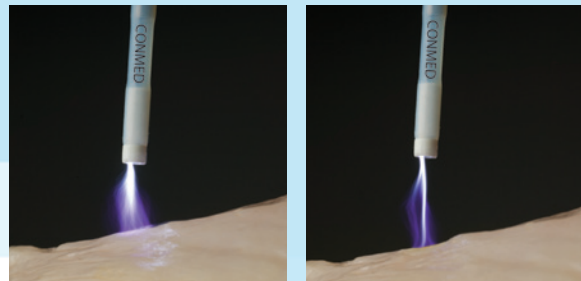

Long, broad and angled beams vary with distance, angle to mucosa and Beamer System settings



Proprietary spiral wire for grabbing flat lesions



Monofilament loop

Beamer Argon Snare Probe

Latex **FREE**

Sterile, Disposable

Catalog Number	Description	Unit
Beamsnare-1	Beamer Snare Probe 2.5MM X 160 (15MM, Firm Spiral Braided Oval Loop) (Bronchoscope)	5/case
Beamsnare-2	Beamer Snare Probe 2.5MM X 230 (15MM, Firm Spiral Braided Oval Loop) (Gastroscope/Colonoscope/Duodenoscopy)	5/case
Beamsnare-3	Beamer Snare Probe 2.5MM X 230 (30MM, Firm Spiral Braided Oval Loop) (Gastroscope/Colonoscope/Duodenoscopy)	5/case
Beamsnare-4	Beamer Snare Probe 2.5MM X 230 (15MM, Oval, Monofilament Firm Loop) (Gastroscope/Colonoscope/Duodenoscopy)	5/case
Beamsnare-5	Beamer Snare Probe 2.5MM X 230 (30MM, Oval, Monofilament Firm Loop) (Gastroscope/Colonoscope/Duodenoscopy)	5/case

CONMED CORPORATION PRODUCT AREAS:

ARTHROSCOPY • ELECTROSURGERY • ENDOSCOPY • ENDOSURGERY • GASTROENTEROLOGY • INTEGRATED SYSTEMS • PATIENT CARE • POWERED INSTRUMENTS • PULMONOLOGY



525 French Road
 Utica, NY 13502
 (315) 797-8375
 Customer Service: 1-800-448-6506
 FAX: 1-800-438-3051
 International Sales: +1(315) 797-8375 • FAX: +1(315) 735-6235
 email: info@mail.conmed.com
 www.conmed.com