

The first cannulated, bioabsorbable,  
independent screw and spiked  
washer for rotator cuff repair.

# BioCuff<sup>®</sup>

**PROVIDING SECURE SOFT TISSUE FIXATION**


- Sutureless design
- Cannulated for simple and quick insertion
- Independent screw and spiked washer
- High pull-out strength
- Manufactured using proprietary self-reinforced polymers




**Linvatec**  
A ConMed Company

# BioCuff<sup>®</sup> C


The BioCuff C implant is a cannulated, bioabsorbable, independent screw and spiked washer designed to provide rigid soft tissue fixation to bone in the repair of rotator cuff tears. The sutureless design eliminates the need for knot tying and provides simple and quick insertion.

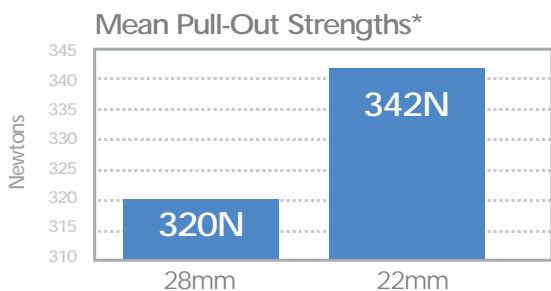


\*Data on file at Linvatec

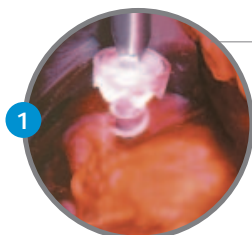

**Hybrid Cortical/Cancellous Thread Design** in multiple lengths provides optimal purchase. (Figure 1)

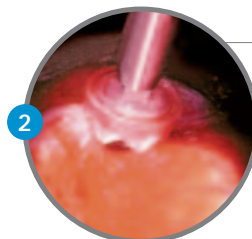

**Independent Spiked Washer** remains stationary as it engages the cuff, eliminating twisting and binding of the tendon as the screw is advanced. (Figure 2)

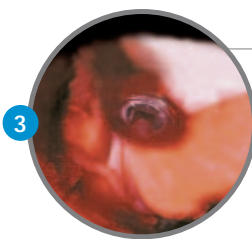

**Sutureless Design** provides rigid and secure fixation without complicated knot tying and suture management. Once fully seated the BioCuff C implant indents the cuff and provides a large area of apposition between the tendon and bone. (Figure 3)



BioCuff C

*Cadaveric testing has shown the BioCuff C implant to have exceptionally high pull-out strength.*



## ORDERING INFORMATION

Ref. No.	Description	Size
236022	BioCuff C (96L/4D) PLA	6.0mm O.D. x 22mm L
236028	BioCuff C (96L/4D) PLA	6.0mm O.D. x 28mm L
236036	BioCuff C (96L/4D) PLA	6.0mm O.D. x 36mm L
Set Consists of:		
SD1002	BioCuff C Cannulated Driver	—
BT1002	BioCuff C Cannulated Soft BoneTap	6.1mm
BT1005	BioCuff C Cannulated Hard Bone Tap	6.3mm
OB1005	BioCuff C Obturator	—
KW1002	BioCuff C K-Wire (6 each)	1.5mm O.D. x 450mm L
CP1002	BioCuff C K-Wire Cap	—
ST1005	BioCuff C Sterilization Tray	—
GU1002	BioCuff C Tissue Grasper	—
CLN3501-04	Cannulated Drill Bit	3.5mm

ARTHROSCOPY • ELECTROSURGERY • ENDOSCOPY • IMAGING • INTEGRATED SYSTEMS • PATIENT CARE • POWERED INSTRUMENTS



11311 Concept Boulevard  
 Largo, FL 33773-4908  
 (727) 392-6464  
 Customer Service: 1-800-237-0169  
 FAX: (727) 399-5256  
 International FAX: +1 (315) 735-6235  
 email: customer\_service@linvatec.com  
 www.linvatec.com

©2003 Linvatec Corporation, a subsidiary of CONMED Corporation, CBR 3024