


STAT 2® I.V. Gravity Flow Controllers

Offering convenience and safety to both patients and healthcare providers

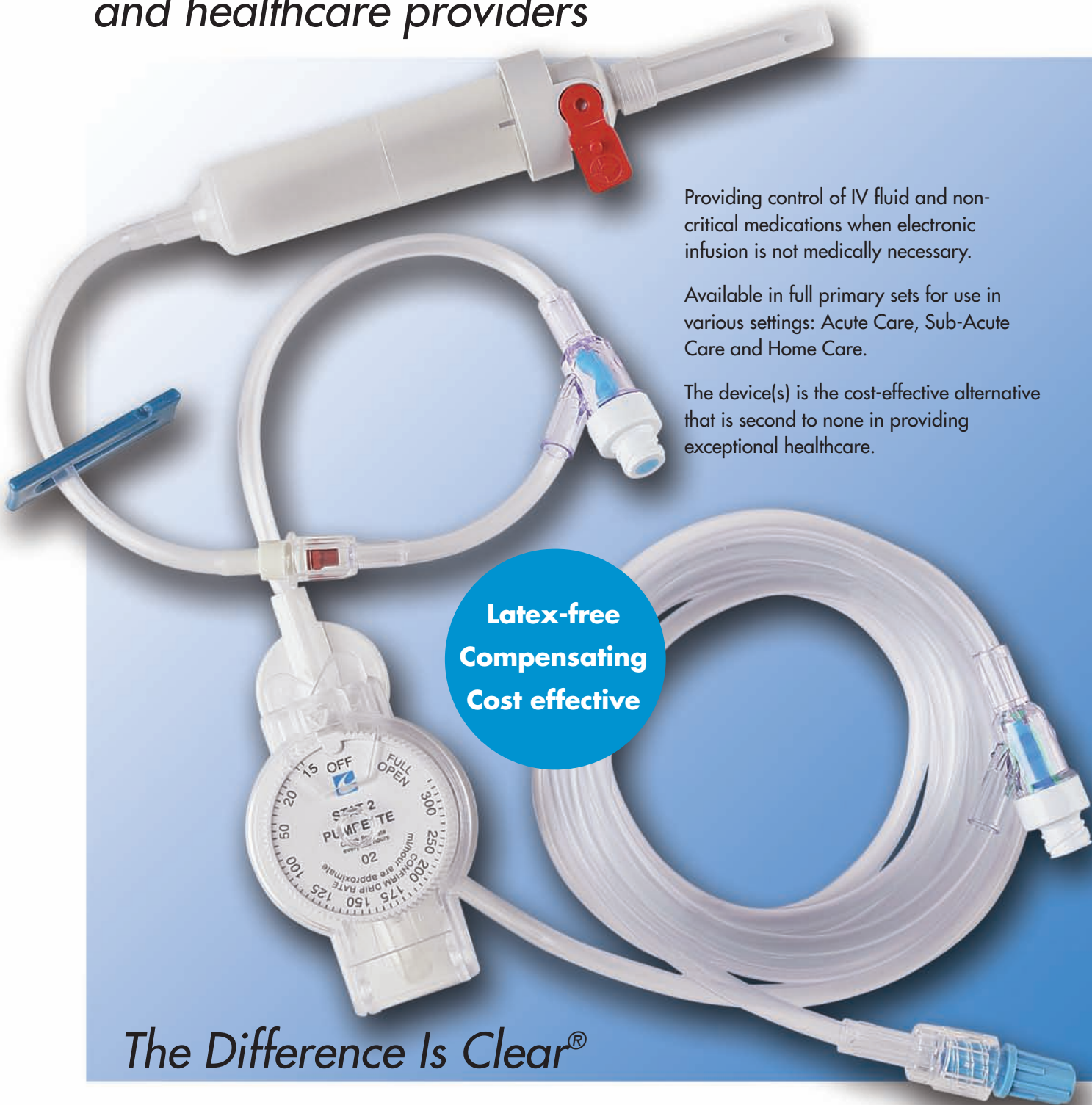

Providing control of IV fluid and non-critical medications when electronic infusion is not medically necessary.

Available in full primary sets for use in various settings: Acute Care, Sub-Acute Care and Home Care.

The device(s) is the cost-effective alternative that is second to none in providing exceptional healthcare.

**Latex-free
Compensating
Cost effective**

The Difference Is Clear®

STAT 2®

I.V. Gravity Flow Control Devices



Designed to control non-critical IV fluid and medications when electronic infusion systems are not necessary.

Offers:

- One handed operation
- Security window with snap lock minimizing accidental manipulation of dial during therapy.
- Ideal for routine hydration and antibiotic therapy

20 and 60-drop primary sets with 15-micron filter within universal drip chamber, back-check valve, a distal luer/slip lock connector and 84" length.

Available with either 2 non-latex split septum "y" sites or 2 non-latex swabbable valve "y" sites


Also available in low profile for extension sets of 12" with 1 non-latex "y" site either split septum or swabbable valve

All IV controllers contain, slide clamp, distal luer/slip lock connection. All primary sets include a back-check valve and universal spike with a 15 micron filter. **All CONMED products are Latex free.**

The Difference Is Clear®

STAT 2[®] Pumpette

Compensating I.V. Gravity Flow Control Devices

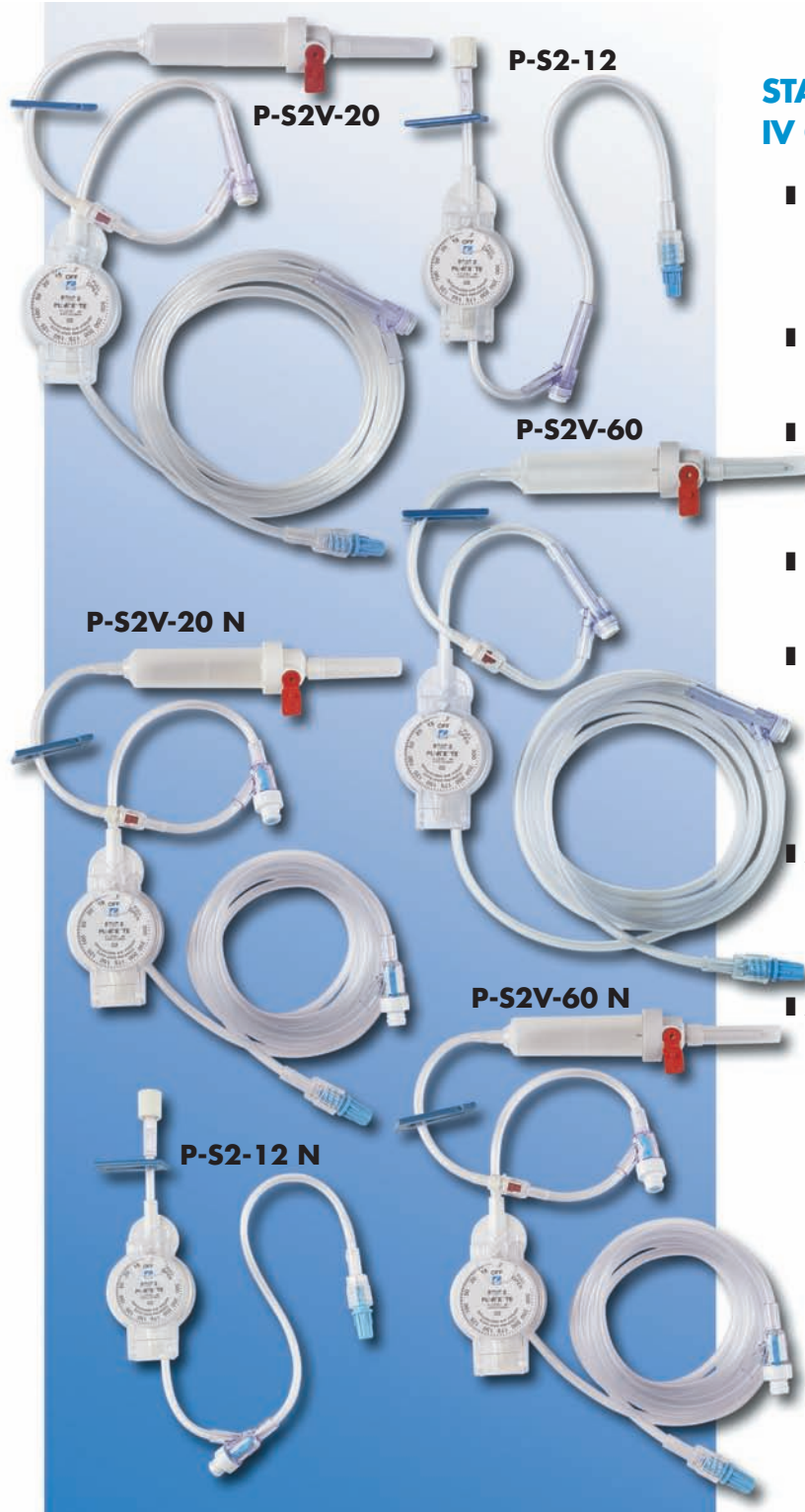

STAT 2[®] Pumpette Compensating IV Controller

■ Compensates for:

- Patients positional changes
- I.V. fluid volumetric changes

■ Controller is at eye level for convenient clinical verification and adjustment

■ Security window minimizes accidental manipulation of dial face during therapy

■ Ideal I.V. delivery system for routine hydration and antibiotic therapies

■ 20 and 60-drop primary sets with 15-micron filter within universal drip chamber, back-check valve, a distal luer/slip lock connector and 84" length.

■ Available with either 2 non-latex split septum "y" sites or 2 non-latex swabbable valve "y" sites

■ Also available in low profile for extension sets of 12" with 1 non-latex "y" site either split septum or swabbable valve

The Difference Is Clear[®]

STAT 2[®] I.V. Gravity Flow Controllers

ORDERING INFORMATION		
STAT 2 [®] CONTROLLERS		
REF. NO.	DESCRIPTION	POUCH/PKG CASE
S2-12	12" length extension set, controller, 1 split-septum "Y" site, distal luer/slip lock connector	1 60
S2V-20	20 drop/ml, 84" primary set controller, 2 split-septum "Y" sites, universal spike, 15 micron-filter, back-check valve, distal luer/slip lock connector	1 48
S2V-60	60 drop/ml, 84" primary set controller, 2 split-septum "Y" sites, universal spike, 15 micron-filter, back-check valve, distal luer/slip lock connector	1 48
S2-12 N	12" length extension set, controller, 1 swabbable valve "Y" site, distal luer/slip lock connector	1 60
S2V-20 N	20 drop/ml, 84" primary set controller, 2 swabbable valve "Y" sites, universal spike, 15 micron-filter, back-check valve, distal luer/slip lock connector	1 48
S2V-60 N	60 drop/ml, 84" primary set controller, 2 swabbable valve "Y" sites, universal spike, 15 micron-filter, back-check valve, distal luer/slip lock connector	1 48
STAT 2 [®] Pumpette Controllers		
P-S2-12	12" length extension set, controller, 1 split-septum "Y" site, distal luer/slip lock connector	1 60
P-S2V-20	20 drop/ml, 84" primary pumpette controller, 2 split-septum "Y" sites, universal spike, 15 micron-filter, back-check valve, distal luer/slip lock connector	1 48
P-S2V-60	60 drop/ml, 84" primary pumpette controller, 2 split-septum "Y" sites, universal spike, 15 micron-filter, back-check valve, distal luer/slip lock connector	1 48
P-S2-12 N	12" length extension set, controller, 1 swabbable valve "Y" site, distal luer/slip lock connector	1 60
P-S2V-20 N	20 drop/ml, 84" primary pumpette controller, 2 swabbable valve "Y" sites, universal spike, 15 micron-filter, back-check valve, distal luer/slip lock connector	1 48
P-S2V-60 N	60 drop/ml, 84" primary pumpette controller, 2 swabbable valve "Y" sites, universal spike, 15 micron-filter, back-check valve, distal luer/slip lock connector	1 48
Accessories		
810-L	Infant size 5.25" x 2.25" arm support w/covers	60 supports, 90 covers
910-L	Pediatric size 7" x 3" arm support w/covers	60 supports, 90 covers
1010-L	Adult size 10.5" x 4" arm support w/covers	80 supports, 80 covers
1011-L	Adult size arm supports ONLY	80 supports
1012-L	Adult size arm support covers	250 covers



525 French Road • Utica, NY 13502
 Customer Service 1-800-448-6506 • FAX 1-800-438-3051
 (315) 797-8375 • FAX (315) 732-7991
 email: info@mail.conmed.com



For International orders or inquiries, please
 contact CONMED International Sales
 (315) 797-8375 • FAX (315) 735-6235
 www.conmed.com