

With both tab and snap options, CONMED can meet all of your needs for the EKG lab. One company with one goal — to satisfy our customer in the most cost-effective way.

## Snaptrace® II

A leader in Resting EKG electrodes, CONMED is proud to introduce the newest addition to the family, the Snaptrace® II Resting Snap electrode. The smaller design is ideal for adults and pediatrics, and offers easy placement for V-leads.

- Silver/Silver Chloride design.
- Solid adhesive hydrogel which provides staying power, but gentle release.
- Radiotranslucent and MRI safe.
- Latex free and PVC free.
- Ideal for short term procedures.

## Fastrace® 4

- Gel 2000 formulation sticks better and stays moist longer — even when the pouch has been opened.
- Silver/Silver Chloride provides stable baselines and excellent electrical characteristics.
- Conductive adhesive gel maintains contact with skin, even on hairy patients.
- Enhanced substrate extends shelf life (36 months).

## HeartBeat

- Gel formulation for improved staying power, but gentle release.
- Smaller size for pediatrics and adults.
- Pliable material is more comfortable and prevents lift off.
- Latex free and PVC free.
- Silver/Silver Chloride provides stable baselines and excellent electrical characteristics.


- Four styles of adapters for HeartBeat and Fastrace® 4.
- Three styles of adapters for Snaptrace® II.
- Adapters compatible with D-Series Individually Shielded Safety Cable Systems:
  - 10 and 14 wire cables available
  - three-year cable warranty.

# RESTING EKG ELECTRODES

# Resting EKG Technical Specifications



Snaptrace® II

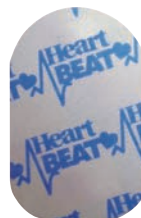

HeartBeat



Fastrace® 4



Parameter	Snaptrace® II	HeartBeat	Fastrace® 4
Reference Number	1921-100	1925-100	1915-100
Quantity	100/pouch; 1000/case	100/pouch; 1000/case	100/pouch; 1000/case
Gel	5% Potassium chloride (KCl) conductive adhesive gel	5% Potassium chloride (KCl) conductive adhesive gel	5% Potassium chloride (KCl) conductive adhesive gel
Gel Surface Area	.9" X 0.65" (22mm X 16mm)	1.25" X 0.812" (32mm X 20mm)	1.375" X 0.875" (35mm X 22mm)
Basepad	Closed cell polyethylene foam	Polyester	Polyester
Basepad Size	0.9" X 1.85" (22mm X 46mm)	1.25" X 0.812" (32mm X 20mm)	1.375 X 0.875" (35mm X 22mm)
Adhesive	Medical grade acrylic adhesive	Conductive adhesive hydrogel	Conductive adhesive hydrogel
Carrier	Polyester	Polyester	Polyester
Shelf life	24 months from date of manufacture	24 months	36 months
Out of pouch shelf life (May vary depending on environmental factors)	8 weeks	30 days	30 days
Special features	Latex free, PVC free, MRI safe Radiotranslucent	Latex free, PVC free Radiotranslucent	Latex free, PVC free Radiotranslucent

HeartBeat and Fastrace® 4 Adapters					Description	
Ref. No.	Fits			Qty/Pkg		
	Banana .156" (4.0mm)	Pin .117" (3.0mm)	Pin .125" (3.2mm)			
ACB	X			10	Alligator clip adapters 6" (15.2cm) long	
ACF		X		10		
ACP			X	10		
HPB-B	X			10	Bi-directional clip adapter 6" (15.2cm) long	
HPF-B		X		10		
HPP-B			X	10		
HPH-Blue	X			10	One piece adapter	
HPH-Red		X	X	10		
<b>Snaptrace® II Adapters</b>					Grabber adapters; 6" (15.2cm) long	
RLB	X			5/10		
RLP			X	5/10		
SNB	X			10		
SNF		X		10	Snap adapters; 6" (15.2cm) long	
SNP			X	10		
<b>Snaptrace® II, HeartBeat and Fastrace® 4 Adapters</b>					Universal adapter	
UMCL	X	X	X	10		

## CONMED CORPORATION PRODUCT AREAS:

ARTHROSCOPY • ELECTROSURGERY • ENDOSCOPY • ENDOSURGERY • GASTROENTEROLOGY • INTEGRATED SYSTEMS • PATIENT CARE • POWERED INSTRUMENTS • PULMONOLOGY



525 French Road  
Utica, NY 13502  
(315) 797-8375

Customer Service: 1-800-448-6506

FAX: 1-800-438-3051

International Sales: +1(315) 797-8375 • FAX: +1(315) 735-6235

email: info@mail.conmed.com

www.conmed.com

© 2009, CONMED Corporation, 04/2009, Control #8328 Rev C

