


P U L S E R E F L E C T A N C E O X I M E T R Y

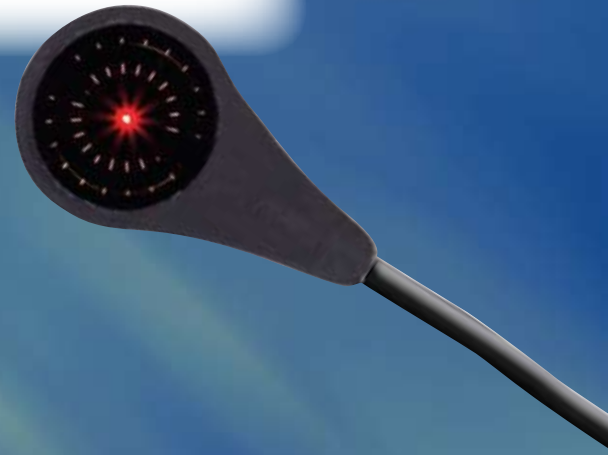

PRO₂[®]

P U L S E R E F L E C T A N C E O X I M E T E R

- **Unique, patented, three wavelength, two detector system:** crosschecks and verifies reading for enhanced performance
- **Flat, low profile sensor:** for use on the forehead or back giving faster response time and greater patient comfort than peripheral sensors
- **Core monitoring:** avoids problems associated with peripheral site monitoring (motion, poor perfusion, etc.)
- **Robust, reusable sensor:** less cost per use compared to disposable sensors
- **For use in adults and pediatrics**

PRO₂[®]

PULSE REFLECTANCE OXIMETER

Unique Features of the PRO₂[®] Pulse Reflectance Oximeter System

Monitor:

- Has a brilliant, easy to read screen
- Turn and press a single knob to navigate all features

Applications:

- Ideal for burn, trauma, hypothermic and hypotensive patients
- Well suited for fragile, friable skin (hydrogel holders, low thermal output, flat sensor)
- When other pulse oximeters have difficulties acquiring accurate values

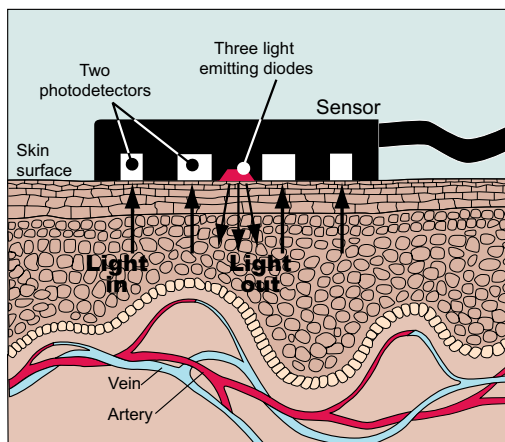


Diagram of PRO₂[®] sensor in operation



- Dimensions** 4.75 inch x 11 inch x 10 inch (12cm x 28cm x 26cm)
- Voltage Input** 100 to 240 VAC at 50 to 60 Hz
- Max AC Power** 40VA consumption
- Battery** Sealed lead-acid for about 2 hours of use

Sensor views (actual size)



Holder views

(shown smaller than actual size)



Top-view PRO-HA/HG



Back-view with Sensor in PRO-HA/HG



Top-view PRO-HS



Back-view PRO-HS

ORDERING INFORMATION	Ref. No.	Product Description	Packaging
	PRO-M	Monitor	1
	PRO-S	Sensor	1
	PRO-HA	Sensor Holder (full-size, acrylic adhesive)	25 ea/box, with 4 boxes/cs
	PRO-HG	Sensor Holder (full-size, hydrogel adhesive)	25 ea/box, with 4 boxes/cs
	PRO-HS	Sensor Holder (small-size, hydrogel adhesive)	25 ea/box, with 4 boxes/cs
	PC-NA	Power Cord	1

Covered by one or more of the following U.S. patents and their international equivalents: 6031603 and 6658277 with other patents pending.

CONMED CORPORATION PRODUCT AREAS:

ARTHROSCOPY • ELECTROSURGERY • ENDOSCOPY • ENDOSURGERY • GASTROENTEROLOGY • INTEGRATED SYSTEMS • PATIENT CARE • POWERED INSTRUMENTS • PULMONOLOGY



525 French Road
 Utica, NY 13502
 (315) 797-8375
 Customer Service: 1-800-448-6506
 FAX: 1-800-438-3051
 International Sales: +1(315) 797-8375 • FAX: +1(315) 735-6235
 email: info@mail.conmed.com
 www.conmed.com

