



ABC[®] DISSECTING ELECTRODES[™]



ABC® Dissecting Electrodes™

ConMed Electrosurgery now provides ABC® technology in contact mode to allow surgeons to dissect while employing the benefits of Argon Beam Coagulation®. Surgeons can also continue to use ABC® in non-contact mode simply by pulling the ABC® Dissecting Electrodes™ away from tissue.

Simple to use, ABC® Dissecting Electrodes™ attach to a ConMed Electrosurgery ABC® handpiece to convert from Argon Beam Coagulation® to Argon Beam Dissect. By forcing argon gas over an electrode in contact or non-contact mode, ABC® Dissecting Electrodes™ may provide the following benefits for surgical procedures:

- **Reduced Smoke from Surgical Field**


Argon gas used with ABC® Dissecting Electrodes™ carries surgical smoke away from the field-of-view.

- **Less Tissue Damage**

Argon gas flowing over the dissection may also help to cool tissue at the perimeter of the intended cut. This may limit thermal spread and may result in less necrosis. The latter is believed to lower the risk of post-operative infections.

- **Clean Cut Hemostasis**

ABC® Dissecting Electrodes™ combine the benefits of argon gas with a tungsten electrode. This combination of electrosurgical modality provides a blend of cut and coagulation effects that may result in a cleaner cut with hemostasis.

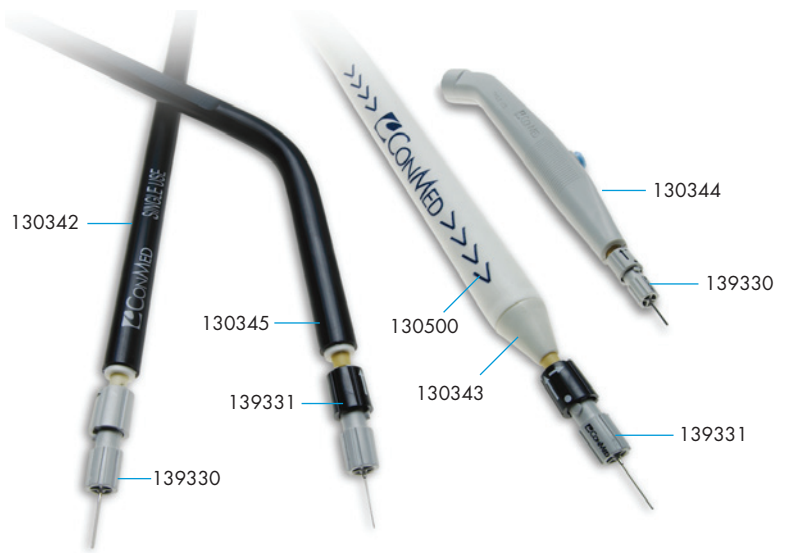

Product	Product Description	Notes	Quantity
139330	ABC® Dissecting Blade Electrode	Converts Argon Beam Coagulation® only to Argon Beam Dissect	20/case
139331	ABC® Dissecting Blunt Needle Electrode	Converts Argon Beam Coagulation® only to Argon Beam Dissect	20/case
130342	Laparoscopic 10 mm ABC® Probe	For use in laparoscopic cases	10/case
130343	ABC® Single Use Nozzle (for use with ABC® Reusable Handpiece, p/n 130500)	For use in open cases. The environmentally friendly choice.	20/case
130344	ABC® Single Function Handpiece	For use in open cases	10/case
130345	ABC® Angled 45° Handpiece	For use in hard-to-reach open cases	10/case

CONMED CORPORATION PRODUCT AREAS:

ARTHROSCOPY • ELECTROSURGERY • ENDOSCOPY • ENDOSURGERY • GASTROENTEROLOGY • INTEGRATED SYSTEMS • PATIENT CARE • POWERED INSTRUMENTS • PULMONOLOGY



525 French Road
 Utica, NY 13502
 (315) 797-8375
 Customer Service: 1-800-448-6506
 FAX: 1-800-438-3051
 International Sales: +1(315) 797-8375 • FAX: +1(315) 735-6235
 email: info@mail.conmed.com
 www.conmed.com

