

# P U L S E R E F L E C T A N C E O X I M E T R Y



## PRO<sub>2</sub><sup>®</sup>

### P U L S E R E F L E C T A N C E O X I M E T E R

- **Unique, patented, three wavelength, two detector system:** crosschecks and verifies reading for improved performance
- **Flat, low profile sensor:** for use on the forehead or back giving faster response time and greater patient comfort
- **Core monitoring:** avoids problems associated with peripheral site monitoring (motion, poor perfusion, etc.)
- **Robust, reusable sensor:** less cost per use compared to disposable sensors
- **For use in adults and pediatrics**

# PRO<sub>2</sub><sup>®</sup>

## PULSE REFLECTANCE OXIMETER

### Unique Features of the PRO<sub>2</sub><sup>®</sup> Pulse Reflectance Oximeter System

#### Monitor:

- Has a brilliant, easy to read screen
- Turn and press a single knob to navigate all features

#### Applications:

- Excellent for stress-test and sleep studies
- Anywhere other pulse oximeters have difficulties acquiring accurate values

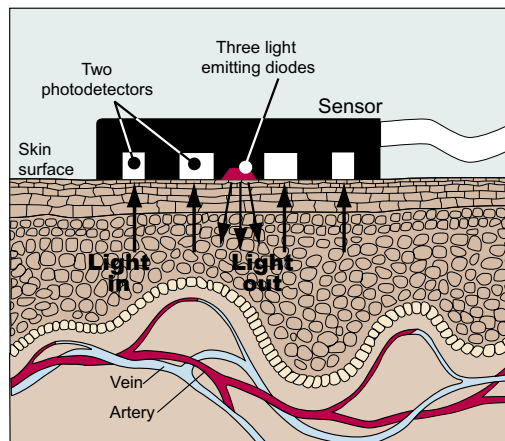


Diagram of PRO<sub>2</sub><sup>®</sup> sensor in operation



PRO-M

**Dimensions** 4.75" x 11" x 10" (12cm x 28cm x 26cm)  
**Voltage Input** 100 to 240 VAC at 50 to 60 Hz  
**Max AC Power** 40VA consumption  
**Battery** Sealed lead-acid for about 2 hours of use



Sensor views (actual size)

Face-view PRO-S



Side-view PRO-S

Sensor holder — Adult  
(shown smaller than actual size)



Top-view PRO-H



Back-view with Sensor PRO-H

ORDERING INFORMATION	Ref. No.	Product Description	Packaging
	PRO-M	Monitor	1
	PRO-S	Sensor	1
	PRO-H	Holder — Adult	25 ea/box, 4 boxes/cs, 100 ea cs
	Call Customer Service	Power cord (country specific)	1

Covered by one or more of the following U.S. patents and their international equivalents: 6031603 and 6658277 with other patents pending.

#### CONMED CORPORATION PRODUCT AREAS:

ARTHROSCOPY • ELECTROSURGERY • ENDOSCOPY • ENDOSURGERY • GASTROENTEROLOGY • INTEGRATED SYSTEMS • PATIENT CARE • POWERED INSTRUMENTS • PULMONOLOGY



525 French Road  
 Utica, NY 13502  
 (315) 797-8375  
 Customer Service: 1-800-448-6506  
 FAX: 1-800-438-3051  
 International Sales: +1(315) 797-8375 • FAX: +1(315) 735-6235  
 email: info@mail.conmed.com  
 www.conmed.com

©2004 ConMed Corporation 9/2005, Control #12451



MBR 9079 Rev 1