


PALADIN™

ROTATOR CUFF ANCHOR

FEEL THE FIXATION STRENGTH OF THE PALADIN ROTATOR CUFF ANCHOR



- 2 Step Insertion Process = Faster Procedure
- Innovative and Aggressive Thread Lead-in on Self Drilling Tap
- High Strength 96L / 4D Self-Reinforced™ Poly Lactic Acid
- 540 N and 440 N Pull-out Strengths
- Low Profile Head of Only 2.75mm
- Proximal and Distal Laser Etch Indicating Eyelet Orientation
- Distal Etch Insertion Depth Guides

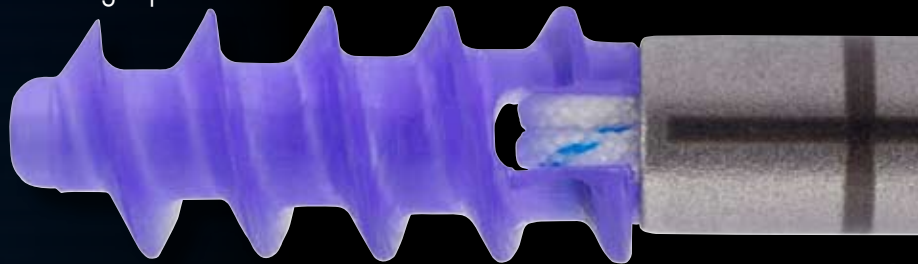

Paladin™ Rotator Cuff Anchor

PALADIN™ ROTATOR CUFF ANCHOR



- The Paladin suture anchors are 5.0mm and 6.5mm diameter screw-in implants manufactured from ConMed Linvatec's patented Self-Reinforced 96L/4D Poly Lactic Acid.
- Each Anchor is pre-loaded on a disposable driver and is pre-threaded with two strands of #2 Hi-Fi® high strength sutures.
- The combination of high pull-out strength and ideal bioabsorbable characteristics in two different sizes will make the Paladin™ anchor an ideal choice for all shoulder rotator cuff procedures.
- Innovative thread lead in design allowing for easy tap engagement

Description

Cat. No.

Paladin 5.0mm Driver with Anchor	C5070H
Paladin 6.5mm Driver with Anchor	C6570H
Paladin 5.0mm Tap	C5075
Paladin 6.5mm Tap	C6575



CONMED CORPORATION PRODUCT AREAS:

ARTHROSCOPY • ELECTROSURGERY • ENDOSCOPY • ENDOSURGERY • GASTROENTEROLOGY • INTEGRATED SYSTEMS • PATIENT CARE • POWERED INSTRUMENTS • PULMONOLOGY



11311 Concept Boulevard
Largo, FL 33773-4908
(727) 392-6464
Customer Service: 1-800-237-0169
FAX: (727) 399-5256
International FAX: +1 (727) 397-4540
email: customer_service@linvatec.com
www.linvatec.com