

UNIVERSAL PLUS[®]

**Electrosurgical and High Volume
Suction/Irrigation Systems for
Surgical Endoscopy**

- Cut
- Cauterize
- Dissect
- Laser
- High Volume Suction
- High Volume Irrigation
- Specimen Retrieval
- Smoke Evacuation

All In One Instrument System



Universal Plus® Straight Handle

A small size pencil grip configuration with color coded E/S buttons or rocker switch activation for precise surgeon control. Familiar trumpet valve S/I buttons are concave/convex and color coded for easy identification. The tethered S/I tubing and E/S cord eliminate entanglements and are directed downward at a 45° angle, to minimize surgeon fatigue.

Electrosurgical Control Options



Hand Controlled Cut and Coag via E/S Buttons with tactile feel and click Hand Controlled S/I



Hand Controlled Cut and Coag Rocker switch Hand Controlled S/I




Foot Controlled Cut and Coag with Universal Cord Hand Controlled S/I

Interchangeable Electrode and Cannula Connection



Universal Plus® Pistol Grip

Familiar pistol grip handle design...for surgeon comfort and control. All activation buttons are color coded. The raised S/I controls offer proportional flow control and the convenient cut/coag buttons feature the desired tactile "click" and S/I tubing and E/S power cord are tethered together to prevent entanglement.

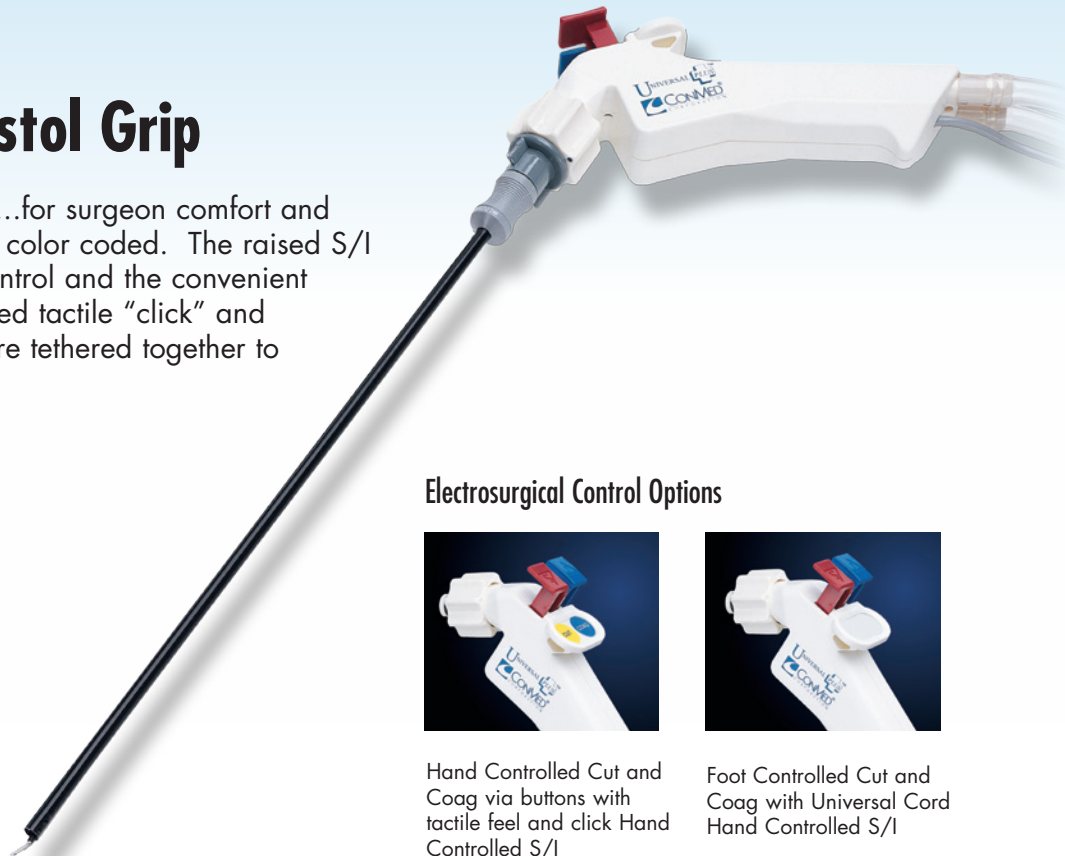

Electrosurgical Control Options



Hand Controlled Cut and Coag via buttons with tactile feel and click Hand Controlled S/I



Foot Controlled Cut and Coag with Universal Cord Hand Controlled S/I

5mm Electrosurgery Electrodes with Suction/Irrigation Optional Stealth ER Eschar Resistant Coating

Available in 27, 32 and 44cm lengths.




Spatula
60-5272-027, 27cm
60-5272-032, 32cm



L-Hook
60-5272-127, 27cm
60-5272-132, 32cm



J-Hook
60-5272-227, 27cm
60-5272-232, 32cm

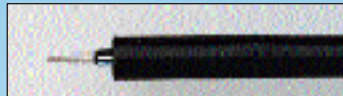

Needle
60-5272-927, 27cm
60-5272-932, 32cm



Spatula
60-5274-027, 27cm
60-5274-032, 32cm
60-5274-044, 44cm



L-Hook
60-5274-127, 27cm
60-5274-132, 32cm
60-5274-144, 44cm




J-Hook
60-5274-227, 27cm
60-5274-232, 32cm
60-5274-244, 44cm




Needle
60-5274-927, 27cm
60-5274-932, 32cm
60-5274-944, 44cm

Extended Insulation at the base of each electrode tip provides greater patient safety, precise surgeon control and reduces possible collateral desiccation. Extendable, lockable tip protector allows coverage of electrode, protecting viscera while suctioning.


Interchangeable Reusable Stainless Steel Suction/Irrigation Cannulas

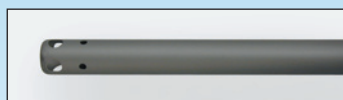


10/10mm Sump Cannula
60-6002-301, Reusable, with side holes, 32cm
60-6002-307, Reusable, with side holes, 15cm

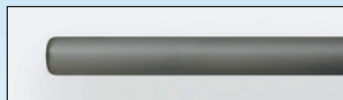

60-6002-306, reusable, no side holes, 32cm
60-6002-305, reusable, no side holes, 15cm

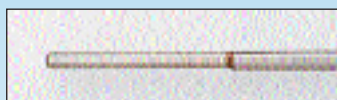


5/3mm Hydrodissector
60-6002-304, Reusable, 32cm

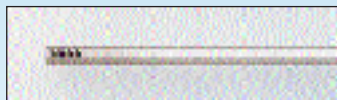

5/5mm Sump Cannula
60-6002-302, reusable, 32cm


Interchangeable Disposable Suction/ Irrigation Cannulas

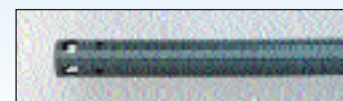


10/10mm Sump Cannula
60-6001-110, Single-Use, 32cm
60-6001-105, Single-Use, no sideholes, 32cm
60-6001-106, Single-Use, no sideholes, 15cm
60-6001-107, Single-Use, with sideholes, 15cm



10/10mm "Soft-Tip" Sump Cannula
60-6000-130, Single-Use, 32cm

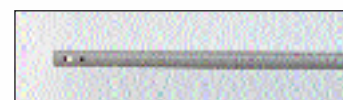

10/5mm Sump Cannula
60-6001-108, Single-Use, 32cm

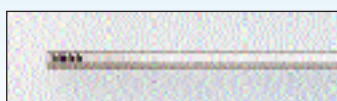


5/5mm Sump Cannula
60-6001-203, Single-Use, 32cm
60-6001-205, Single-Use, 44cm




5/5mm Laser Access Cannula
60-6001-360, Single-Use, 32cm
60-6001-361, Single-Use, 44cm

All with Unique Extendable Lockable Tip Protector Feature

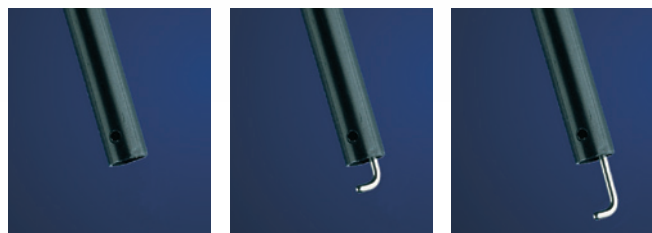

This feature allows varied tip exposure during electrosurgery. During Suction/Irrigation, it protects viscera while suctioning and directs flow while irrigating.

"Dual-Spike" and Standard Suction Connections

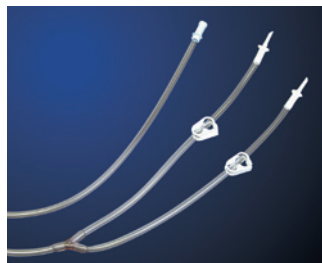

All CONMED Universal Plus® Control Handles feature the standard push-on connector for suction canister attachment, and "dual-spike" irrigation line for quick bag changes or dual pump set-ups.

Ordering Information

Product Number	Description	Packaging	Product Number	Description	Packaging
Control Handles			Electrodes with S/I Channel, Extendable, Lockable Tip Protector and Extended Insulation Optional Stealth ER Eschar Resistant Coating.		
60-6010-000	Pistol Grip, Foot Control E/S, Hand Control S/I, w/10' HiFlo tubing, Irrigation bag spike and universal cord	1/pkg 10/cs	60-5274-132	"L" Hook w/Stealth ER, 32cm, 5mm	1/pkg 5/cs
60-6010-002	Pistol Grip, Hand Control E/S and S/I, w/10' HiFlo tubing, Irrigation bag spike	1/pkg 10/cs	60-5274-144	"L" Hook w/Stealth ER, 44cm, 5mm	1/pkg 5/cs
60-6010-020	Pistol Grip, Foot Control E/S, Hand Control S/I, w/10' HiFlo tubing "Y" Connector for mechanical pumps and universal cord	1/pkg 10/cs	60-5272-227	"J" Hook, 27cm, 5mm	1/pkg 5/cs
60-6010-022	Pistol Grip, Hand Control E/S and S/I, w/10' HiFlo tubing Irrigation "Y" Connector for mechanical pumps	1/pkg 10/cs	60-5272-232	"J" Hook, 32cm, 5mm	1/pkg 5/cs
60-6010-001	Straight Handle, Foot Control E/S and S/I, w/10' HiFlo tubing, Irrig. bag spike and universal cord	1/pkg 10/cs	60-5274-227	"J" Hook w/Stealth ER, 27cm, 5mm	1/pkg 5/cs
60-6010-003	Straight Handle, Hand Control E/S (Buttons) and S/I, w/10' HiFlo tubing and Irrig. bag spike	1/pkg 10/cs	60-5274-232	"J" Hook w/Stealth ER, 32cm, 5mm	1/pkg 5/cs
60-6010-005	Straight Handle, Hand Control E/S (Rocker switch) and S/I, w/10' HiFlo tubing and Irrig. bag spike	1/pkg 10/cs	60-5274-244	"J" Hook w/Stealth ER, 44cm, 5mm	1/pkg 5/cs
60-6020-100	Pistol Grip, Hand Control S/I only w/10' HiFlo tubing, bag spike	1/pkg 10/cs	60-5272-927	Needle, 27cm, 5mm	1/pkg 5/cs
60-6020-115	Straight Handle, Hand Control S/I only w/10' HiFlo tubing, bag spike	1/pkg 10/cs	60-5272-932	Needle, 32cm, 5mm	1/pkg 5/cs
Electrodes with S/I Channel, Extendable, Lockable Tip Protector and Extended Insulation Optional Stealth ER Eschar Resistant Coating			Cannulas-Disposable		
60-5272-027	Spatula, 27cm, 5mm	1/pkg 5/cs	60-6001-203	5mm, Sump, 32cm, Single-Use	1/pkg 5/cs
60-5272-032	Spatula, 32cm, 5mm	1/pkg 5/cs	60-6001-205	5mm, Sump, 44cm, Single-Use	1/pkg 5/cs
60-5274-027	Spatula w/Stealth ER, 27cm, 5mm	1/pkg 5/cs	60-6001-108	10/5mm Sump, 32cm, Single-Use	1/pkg 5/cs
60-5274-032	Spatula w/Stealth ER, 32cm, 5mm	1/pkg 5/cs	60-6000-130	10mm "Soft Tip" sump, 32cm, Single-Use	1/pkg 5/cs
60-5274-044	Spatula w/Stealth ER, 44cm, 5mm	1/pkg 5/cs	60-6001-110	10mm Sump, 32cm, Single-Use	1/pkg 5/cs
60-5272-127	"L" Hook, 27cm, 5mm	1/pkg 5/cs	60-6001-105	10mm Sump, no side holes, 32cm, Single-Use	1/pkg 5/cs
60-5272-132	"L" Hook, 32cm, 5mm	1/pkg 5/cs	60-6001-106	10mm Sump, no side holes, 15cm, Single-Use	1/pkg 5/cs
60-5274-127	"L" Hook w/Stealth ER, 27cm, 5mm	1/pkg 5/cs	60-6001-107	10mm Sump, w/side holes, 32cm, Single-Use	1/pkg 5/cs
			60-6001-360	5mm Laser Cannula, 32cm, Single-Use	1/pkg 5/cs
			60-6001-361	5mm Laser Cannula, 44cm, Single-Use	1/pkg 5/cs
			Cannulas-Reusable Stainless Steel		
			60-6002-302	5mm, 32cm, reusable	1/pkg 1/cs
			60-6002-304	5/3mm Hydrodissector, 32cm, reusable	1/pkg 1/cs
			60-6002-305	10mm Sump, no side holes 15cm, reusable	1/pkg 1/cs
			60-6002-306	10mm Sump, no side holes 32cm, reusable	1/pkg 1/cs
			60-6002-307	10mm Sump, w/side holes 15cm, reusable	1/pkg 1/cs
			60-6002-301	10mm Sump, w/side holes 32cm, reusable	1/pkg 1/cs
			60-5168-005	Universal 10' Monopolar E/S cord, Single-Use	1/pkg 50/cs
			60-5177-003	"Y" Connector set for mechanical irrigation pumps	1/pkg 10/cs
			87101CS	"Y" Set w/dual spike	1/pkg 10/cs

Note: Custom Procedure Packs Available. All products are individually packaged and Sterile for Single Patient Use except the Reusable Cannulas.

CONMED CORPORATION PRODUCT AREAS:

ARTHROSCOPY • ELECTROSURGERY • ENDOSCOPY • ENDOSURGERY • GASTROENTEROLOGY • INTEGRATED SYSTEMS • PATIENT CARE • POWERED INSTRUMENTS • PULMONOLOGY



525 French Road
 Utica, NY 13502
 (315) 797-8375
 Customer Service: 1-800-448-6506
 FAX: 1-800-438-3051
 International Sales: +1(315) 797-8375 • FAX: +1(315) 735-6235
 email: info@mail.conmed.com
 www.conmed.com