



NEW

Cannula Portfolio


Part of CONMED's Shoulder Restoration System™ For arthroscopic shoulder procedures, CONMED has a diverse portfolio of unique cannulas in multiple sizes, materials, and thread designs to meet a wide range of needs and applications.



Y-KNOT™ FLEX
ALL-SUTURE ANCHOR



POPLOK®
KNOTLESS ANCHOR

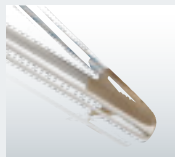

SPECTRUM™ AUTOPASS™
SUTURE PASSER



To learn more about these and other innovative products, call **800-237-0169** or visit **CONMED.com**.

SHOULDER
RESTORATION
SYSTEM™

ARC™

Adjustable Retractor Cannula

A true innovation, the unique design of CONMED's ARC Adjustable Retractor Cannula helps surgeons see more clearly, operate with greater range of motion and minimize cannula pull-out.



Better Visualization

The proprietary locking nut and distal wing design of this cannula retracts the deltoid fascia to create a large space for better visualization.

More Room in the Subacromial Space


This muscle retraction feature also creates more room in the subacromial space – making it easier to use and manipulate elongated instruments like suture passers.

Secure, Customizable Fit

The unique locking nut and wing design allows surgeons to optimize the fit to each individual patient while securely locking the cannula in place.

Unique Nut and Wing Design

Minimizes cannula pull-out and avoid a loss of distension.

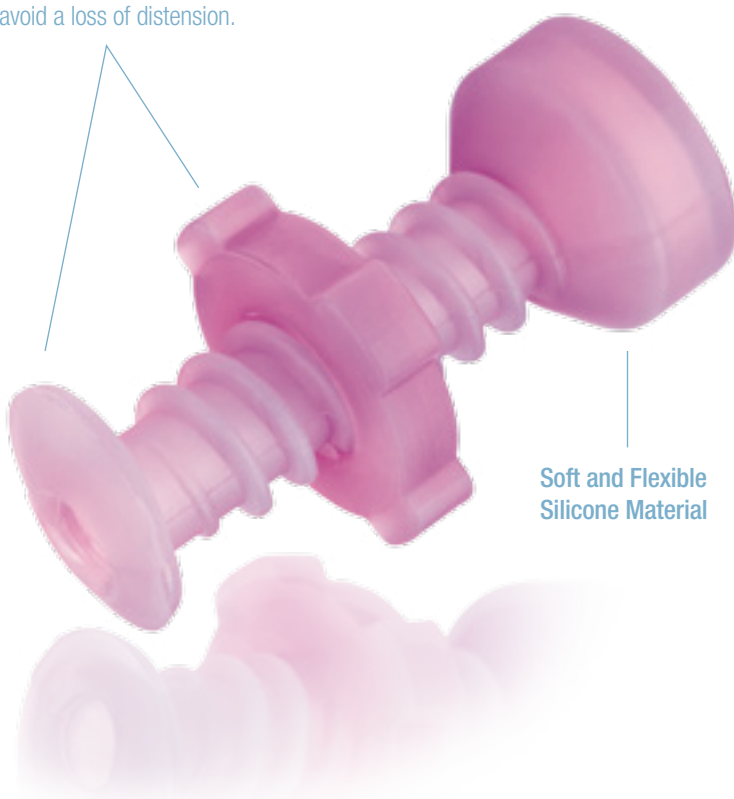

Soft and Flexible
Silicone Material

Full Range of Motion

In addition to helping protect tissue and cartilage, the soft and flexible silicone material allows surgeons to more easily handle and manipulate instruments. This helps place instruments at angles and locations that would be challenging with a rigid cannula.

SHOULDER
RESTORATION
SYSTEM™



20° Acute Angle Wing


SHOULDER
RESTORATION
SYSTEM™



Rounded Tip

Dry-Doc® Cannula

The Dry-Doc® cannula is an all-in-one smooth and ribbed cannula and provides a robust hold in skin, muscle and capsule. It has double tri-cut seals, a click-lock ratcheting stopcock for precise fluid control, and a rounded distal tip to prevent damage to cartilage, soft tissue and suture. The Dry-Doc cannula is translucent for unobstructed visualization and is available in multiple sizes and lengths.

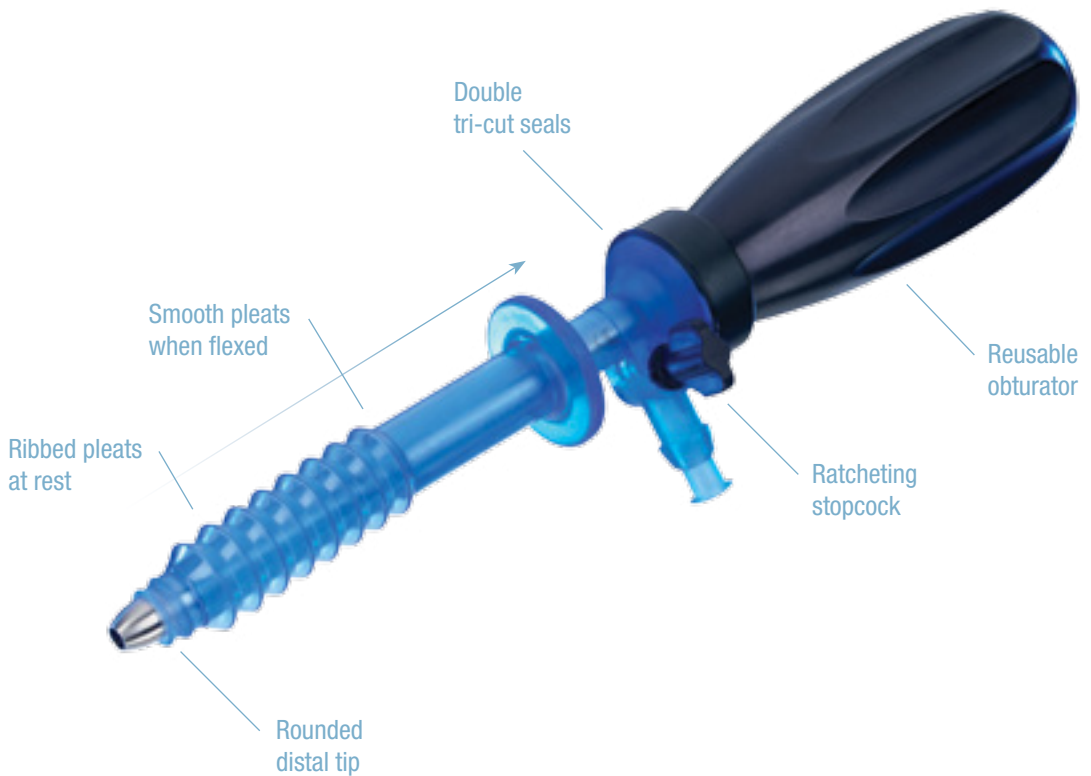



Hex Flex[®] Cannula

Pliable, bendable, and ultra-durable, the Hex Flex[®] cannula features two distal spiral threads for a sturdy yet atraumatic hold in capsule. The smooth shaft allows easy manipulation, while the hex cap interface provides easier insertion into the joint capsule. Double tri-cut seals prevent leaks and the click-lock ratcheting stopcock permits precise fluid control. The Hex Flex is translucent for unobstructed visualization and is available in multiple sizes.

SHOULDER
RESTORATION
SYSTEM™




**Hex Interface
with Obturator**

SHOULDER RESTORATION SYSTEM™



Double Seal

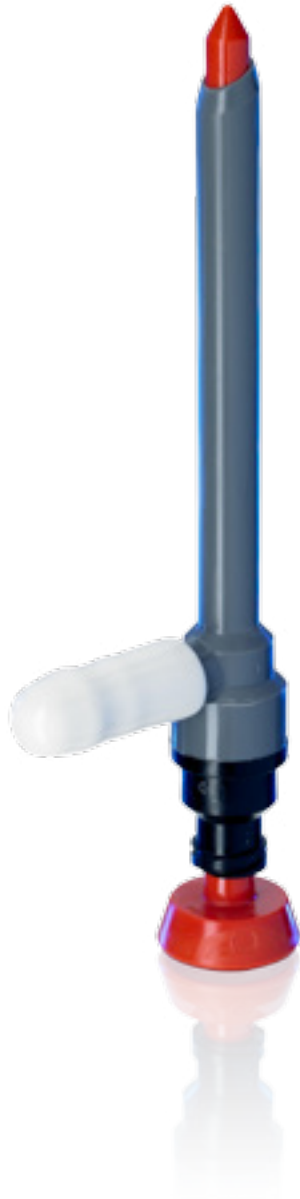

Universal Cannula

Versatile for a wide range of needs, CONMED's Universal Cannula Set features a double seal to help prevent backwash of fluids. The connector allows positive connection to some CONMED blades and burs and the material is non-conductive for electro-surgical procedures. Sharp and blunt trocars are included in the set.



Clear Flexible Cannula

This cannula's flexibility and transparent design enable surgeons to use larger instruments and curved blades, as well as make it easier to tie knots. Available in threaded or smooth designs.



Rigid Cannulas

A versatile, all-purpose solution, CONMED's Rigid Disposable cannulas feature a translucent body design that allows visualization of passing instrumentation or suture. The dual-silicone seal system and integrated flow control valve also help prevent fluid leakage. For convenience, both luer and 1/4" tube connection are on a single fluid port and a disposable conical obturator is included as well. Available in threaded or smooth designs.

SHOULDER
RESTORATION
SYSTEM™



Fully Threaded



525 French Road, Utica, NY 13502

SHOULDER RESTORATION SYSTEM™

Ordering Information

	Description	Size	Catalog Number
NEW	ARC™ Cannula		
	ARC Cannula, green	6 x 65mm	C06-65
	ARC Cannula with outflow, green	6 x 85mm	C06-85
	ARC Cannula, purple	8 x 65mm	C08-85
	ARC Cannula with outflow, purple	8 x 85mm	C08-85
	14cm Forceps, sterile	14cm	SU-9-0183-14
	Dry-Doc® Cannula		
	Dry-Doc Cannula, translucent blue	7 x 85mm	C7360
	Reusable Obturator, anodized blue (for use with C7360)	7 x 85mm	C7385
	Dry-Doc Cannula, translucent blue w/ black cap	7 x 95mm	C7369
	Reusable Obturator, anodized black (for use with C7369)	7 x 95mm	C7388
	Hex Flex® Cannula		
	Hex Flex Cannula, translucent green	5 x 85mm	C7450
	Reusable Obturator, anodized green (for use with C7450)	5 x 85mm	C7380
	Hex Flex Cannula, translucent gray	6 x 85mm	C7455
	Reusable Obturator, anodized gray (for use with C7455)	6 x 85mm	C7382
	Universal Cannula		
	Universal Cannula set with fenestrations	5.5 x 70mm	9704
	Universal cannula set without fenestrations	5.5 x 70mm	9718
	Clear Flexible Cannula		
	Clear flexible cannula set, threaded	6.5 x 73mm	C7312
	Clear flexible cannula set, smooth	6.5 x 73mm	C7313
	Reusable cannulated obturator		9828
	Legacy Cannula		
	Legacy threaded cannula set	6 x 50mm	C7322
		8.4 x 90mm	C7372

NOTE: More sizes and lengths available at CONMED.com

Local 727-392-6464
Toll Free 1-800-237-0169

CONMED.com
customer_service@conmed.com