

Stealth

Diamond-Like Carbon Coated Shaver Blades

CONMED's top-of-the-line performance blades, Stealth arthroscopic shaver blades feature a hard, gall-resistant Diamond-Like Carbon coating to reduce metal particulate generation,^{1,2} lubrication on the inner tube to reduce metal-on-metal contact, and our patented XtraSharp™ technology for optimum cutting performance.



Features and Benefits

Reduced Metal Particulate Generation

Diamond-Like Carbon – coating on distal cutting tip has a lower coefficient of friction than standard steel and has shown reduced metal particulate generation when compared to leading competitors.^{1,2}

Lubrication – (full length shrink coverage) on the inner tube provides a barrier designed to prevent metal-on-metal contact – further helping reduce metal particulate generation within the blade shaft.

Cleaner, More Efficient Cut

Two-piece inner cutting tip – design provides a 22% tighter clearance between the inner and outer cutting window allowing for a cleaner, more efficient cut.³

Optimum Cutting Performance


Patented XtraSharp™ – electrochemical grinding method provides a sharp cutting edge for ideal cutting performance.

Stealth's Rigid TrimRite® Cutting Tip – is harder than standard Nitronic® 60 grade steel, which is designed to improve cutting efficiency and reduce particulate generation.³

To learn more about this and other innovative products, call **1-800-237-0169** or visit **CONMED.com**.

The design features of the Stealth DLC Coated Shaver Blades help to reduce metal particulate generation against a leading competitor – as shown in the data comparison chart below.

Metal Particulate Comparison Chart*

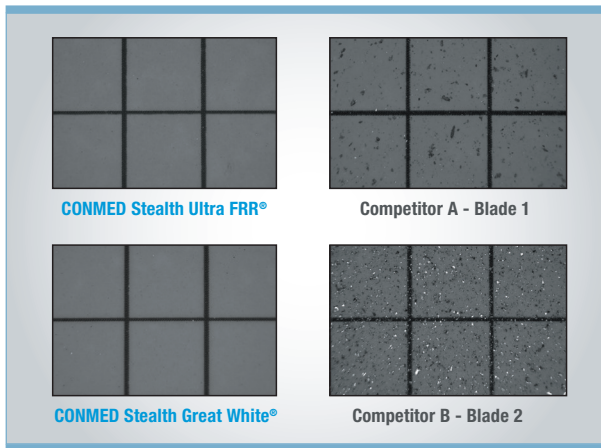

*TR17-00073, TR17-00165. Competitor A - Blade 1: Smith & Nephew, Dyonics Bonecutter. Competitor A - Blade 2: Smith & Nephew, Incisor.

UltraFRR®
Stealth DLC Coated Blade

Great White®
Stealth DLC Coated Blade



Ordering Information

Description	Size	Length	Qty.	Catalog No.
Stealth Diamond-Like Carbon (DLC) Coated Blades				
Stealth Blade, Great White®	4.2mm	13cm	1	S9299A
Stealth Blade, UltraFRR®	4.2mm	13cm	1	SDPS-C010

CONMED Corporation
525 French Road, Utica, NY 13502

Local 727-392-6464
Toll Free 1-800-237-0169

CONMED.com
customer_service@conmed.com