



# Shoulder Cannula Portfolio

For arthroscopic shoulder procedures, CONMED has a diverse portfolio of unique cannulas in multiple sizes, materials, and thread designs to meet a wide range of needs and applications.



Learn more about HydroSeal™  
& other cannula options call  
**1-866-4CONMED** or visit  
[www.CONMED.com](http://www.CONMED.com).



Imagine what we could do together.

# HydroSeal™ Silicone Cannula

## FLEXIBLE SILICONE CANNULAS

Experience seamless access with CONMED's flexible silicone cannula - enhancing maneuverability, safeguarding soft tissues, and maintaining a clear surgical field.

### Fluid Control

The double dam helps to create a more effective seal, reducing fluid leakage during the procedure. By controlling fluid flow, the double dam design reduces fluid leakage.

### Adjustable Stability

Ribs along the axis provide additional stability by preventing the cannula from sliding in and out of the incision. The ribs allow the cannula length to be adjusted, accommodating variations in patient anatomy and surgical needs.



## Optimal Fit

The removable spacer optimizes the fit within the portal, enhancing stability and ease of use during the procedure. In addition, it provides rigidity to easily shuttle instruments through the cannula.

## Low Profile Distal Flange

The distal flange is designed to sit flush against soft tissue, so that it will not interfere with the workspace anatomy and surgical needs.



Imagine what we could do together.

## Dry-Doc<sup>®</sup> Silicone Cannula

The Dry-Doc<sup>®</sup> cannula is an all-in-one smooth and ribbed cannula that provides a robust hold in skin, muscle, and capsule. It has double tri-cut seals, a click-lock ratcheting stopcock for precise fluid control, and a rounded distal tip to prevent damage to cartilage, soft tissue, and suture. The Dry-Doc<sup>®</sup> cannula is translucent for unobstructed visualization and is available in multiple sizes and lengths.



Hex Interface  
with Obturator

## Hex Flex<sup>®</sup> Cannula


Pliable, bendable, and ultra-durable, the Hex Flex<sup>®</sup> cannula features two distal spiral threads for a sturdy yet atraumatic hold in capsule. The smooth shaft allows easy manipulation, while the hex cap interface provides easier insertion into the joint capsule. Double tri-cut seals prevent leaks and the click-lock ratcheting stopcock permits precise fluid control. The Hex Flex<sup>®</sup> is translucent for unobstructed visualization and is available in multiple sizes.

## Universal Cannula

Versatile for a wide range of needs, CONMED's Universal Cannula Set features a double seal to help prevent backwash of fluids. The connector allows positive connection to some CONMED blades and burrs, and the material is non-conductive for electrosurgical procedures. Sharp and blunt trocars are included in the set.



Double Seal

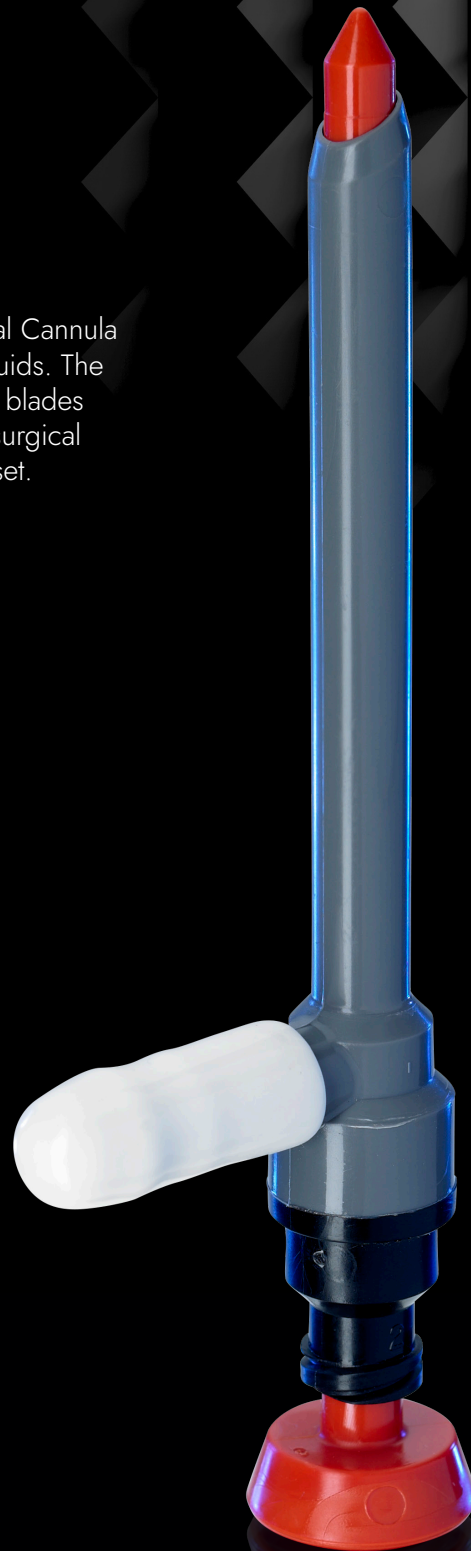

## Clear Flexible Cannula

This cannula's flexibility and transparent design enable surgeons to use larger instruments and curved blades, as well as make it easier to tie knots. Available in threaded or smooth designs.



Fully Threaded



## Rigid Cannulas

A versatile, all-purpose solution, CONMED's Rigid Disposable cannulas feature a translucent body design that allows visualization of passing instrumentation or suture. The dual-silicone seal system and integrated flow control valve also help prevent fluid leakage. For convenience, both luer and 1/4" tube connection are on a single fluid port, and a disposable conical obturator is also included. Available in threaded or smooth designs.



# Cannula Portfolio

## Ordering Information

Description	Size	Catalog Number
<b>Dry-Doc® Cannula</b>		
Dry-Doc® Cannula, Translucent Green, Packaged Sterile, with a Disposable Obturator	5 x 85mm	C7350
Dry-Doc® Cannula, Translucent Gray, Packaged Sterile, with a Disposable Obturator	6 x 85mm	C7359
Dry-Doc® Cannula, Translucent Blue, Packaged Sterile, with a Disposable Obturator	7 x 85mm	C7360
Dry-Doc® Cannula, Translucent Red, Packaged Sterile, with a Disposable Obturator	8 x 75mm	C7367
Dry-Doc® Cannula, Translucent Yellow, Packaged Sterile, with a Disposable Obturator	8 x 85mm	C7368
Dry-Doc® Cannula, Translucent Blue w/ Black cap, Packaged Sterile, with a Disposable Obturator	7 x 95mm	C7369
Reusable Obturator, anodized blue (for use with C7360)	7 x 85mm	C7385
Reusable Obturator, anodized black (for use with C7369)	7 x 95mm	C7388
<b>Hex Flex® Cannula</b>		
Reusable Obturator, anodized green (for use with C7450)	5 x 85mm	C7380
Hex Flex® Cannula, Translucent Green, Packaged Sterile, with Disposable Obturator	5 x 85mm	C7450
Hex Flex® Cannula, Translucent Gray, Packaged Sterile, with Disposable Obturator	6 x 85mm	C7455
Hex Flex® Cannula, Translucent Blue, Packaged Sterile, with Disposable Obturator	7 x 85mm	C7465
Hex Flex® Cannula, Translucent Yellow, Packaged Sterile, with Disposable Obturator	8 x 85mm	C7480
Reusable Obturator, anodized gray (for use with C7455)	6 x 85mm	C7382
<b>Universal Cannula</b>		
Universal Cannula set with fenestrations	5.5 x 70mm	9704
Universal cannula set without fenestrations	5.5 x 70mm	9718
<b>Clear Flexible Cannula</b>		
Clear flexible cannula set, threaded	6.5 x 73mm	C7312
Clear flexible cannula set, smooth	6.5 x 73mm	C7313
Reusable cannulated obturator		9828



# Cannula Portfolio

## Ordering Information

Description	Size	Catalog Number
<b>Rigid Cannula</b>		
Threaded Disposable Cannula with Disposable Obturator	6 x 50mm	C7322
Smooth Disposable Cannula with Disposable Obturator	6 x 50mm	C7324
Threaded Disposable Cannula with Disposable Obturator	6 x 75mm	C7332
Smooth Disposable Cannula with Disposable Obturator	6 x 75mm	C7334
Threaded Disposable Cannula with Disposable Obturator	6 x 90mm	C7342
Smooth Disposable Cannula with Disposable Obturator	8.4 x 50mm	C7352
Threaded Disposable Cannula with Disposable Obturator	8.4 x 50mm	C7354
Smooth Disposable Cannula with Disposable Obturator	8.4 x 75mm	C7362
Threaded Disposable Cannula with Disposable Obturator	8.4 x 75mm	C7364
Threaded Disposable Cannula with Disposable Obturator	8.4 x 90mm	C7372
Smooth Disposable Cannula with Disposable Obturator	8.4 x 90mm	C7374
<b>HydroSeal™ Silicone Cannula</b>		
HydroSeal™ Silicone Cannula	8x20mm	HS-0820
HydroSeal™ Silicone Cannula	8x30mm	HS-0830
HydroSeal™ Silicone Cannula	8x40mm	HS-0840
HydroSeal™ Silicone Cannula	8x50mm	HS-0850
HydroSeal™ Silicone Cannula	10x30mm	HS-1030
HydroSeal™ Silicone Cannula	10x40mm	HS-1040
HydroSeal™ Silicone Cannula	10x50mm	HS-1050
HydroSeal™ Silicone Cannula	12x40mm	HS-1240

CONMED Corporation | 11311 Concept Blvd., Largo, FL 33773 | [www.CONMED.com](http://www.CONMED.com)



Service & Support