

Eliminator® PET Balloon Dilators

CONMED's Eliminator® PET Balloon Dilators offer a simple solution to dilation even in the tightest of strictures. Constructed of patented PET material, these dilators deploy a multi-stage inflation technique that eliminates the need for more than one balloon.












To learn more about these and other innovative products, call 1-866-4CONMED or visit www.CONMED.com

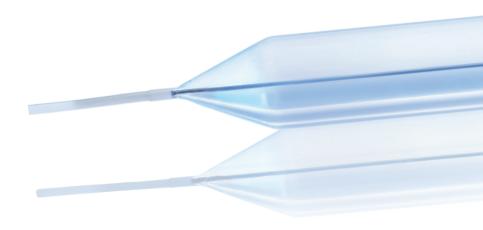

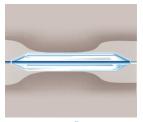

ELIMINATOR® PET BALLOON DILATOR

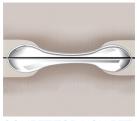

COMPETITOR NON PET BALLOON DILATOR

Pressure Dilates Strictures, Not Size

Balloons dilate strictures with the radial force generated by inflating them to predetermined pressures. The Eliminator® balloons patented PET material disperses pressure along the full length of the balloon ensuring repeatable, reliable balloon sizing. When resistance is met with non PET balloons, the balloon may expand into the area surrounding the stricture. When this occurs, the expected balloon diameter may not be achieved when inflated to recommended levels.


In contrast, Eliminator PET Balloon Dilators are true to size— every time

Balloon Size	Stage One	Stage Two	Stage Three*
6mm /18F	60psi/4.1atm	90psi/6.1atm	120psi/8.2atm
8mm /24F	60psi/4.1atm	90psi/6.1atm	120psi/8.2atm
10mm /30F	45psi/3.1atm	68psi/4.6atm	90psi/6.1atm
12mm /36F	45psi/3.1atm	68psi/4.6atm	90psi/6.1atm
15mm /45F	30psi/2.0atm	45psi/3.1atm	60psi/4.1atm
18mm /54F	25psi/1.7atm	38psi/2.6atm	50psi/3.4atm
20mm /60F	25psi/1.7atm	38psi/2.6atm	50psi/3.4atm

The Eliminator® PET Balloon Advantage

- Accurate Sizing at Maximum Pressure
- Uniform Expansion/No Waisting
- Valued Performance

^{*}Device should not exceed maximum inflation pressure defined by stage three.

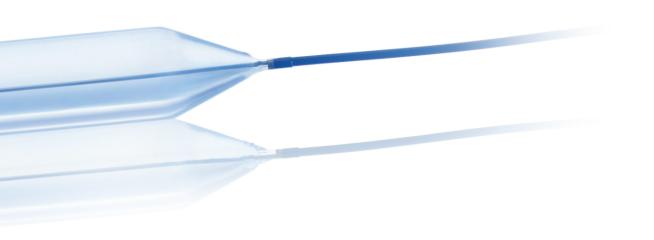

CONMED's Rule of One

The CONMED Eliminator® PET Balloon Dilators combine the radial strength you need to dilate even the tightest strictures with an innovative and cost-effective alternative. The Multi-Stage Inflation Technique is a simple method for obtaining a complete and reliable dilation with a single Eliminator PET Balloon Dilator.

A Cost-Effective Approach


CONMED aims to provide high quality products at an affordable price. In alignment with this goal, Eliminator® PET Balloon Dilator is one of the more cost-effective multi-stage balloon dilation systems in the GI space. Based on 2020 market data, the average price for Eliminator® is approximately 59% less than the average industry price for multi-stage balloon dilators. Enjoy an appreciable cost savings without compromising on quality or customer service.





BALLOON INFLATION DEVICE

Average Unit Sale Price

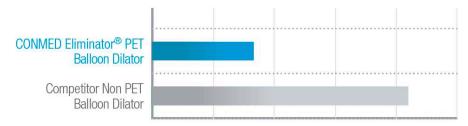





Description	Balloon Diameter		Catheter Length	Unit	Catalog Numbe
Eliminator® PET Biliary Balloon Dilators					
Esophageal Balloon Dilator	6mm	8cm	180cm	1/box	000340
Esophageal Balloon Dilator	8mm	8cm	180cm	1/box	000341
Esophageal Balloon Dilator	10mm	8cm	180cm	1/box	000342
Esophageal Balloon Dilator	12mm	8cm	180cm	1/box	000343
Esophageal Balloon Dilator	15mm	8cm	180cm	1/box	000344
Esophageal Balloon Dilator	18mm	8cm	180cm	1/box	000345
Esophageal Balloon Dilator	20mm	8cm	180cm	1/box	000380
Esophageal Advantage Pack (5ea. of No.000344)	15mm	8cm	180cm	5/box	000735
Esophageal Advantage Pack (5ea. of No.000345)	18mm	8cm	180cm	5/box	000736
Esophageal Advantage Pack (5ea. of No.000380)	20mm	8cm	180cm	5/box	000379
Pyloric/Colonic Balloon Dilator	6mm	4cm	230cm	1/box	000844
Pyloric/Colonic Balloon Dilator	8mm	4cm	230cm	1/box	000845
Pyloric/Colonic Balloon Dilator	10mm	4cm	230cm	1/box	000846
Pyloric/Colonic Balloon Dilator	12mm	4cm	230cm	1/box	000847
Pyloric/Colonic Balloon Dilator	15mm	4cm	230cm	1/box	000848
Pyloric/Colonic Balloon Dilator	18mm	4cm	230cm	1/box	000849
Pyloric/Colonic Balloon Dilator	20mm	4cm	230cm	1/box	000861
Pyloric/Colonic Advantage Pack (5ea. of No.000848)	15mm	4cm	230cm	5/box	000858
Pyloric/Colonic Advantage Pack (5ea. of No.000849)	18mm	4cm	230cm	5/box	000859
Pyloric/Colonic Advantage Pack (5ea. of No.000861)	20mm	4cm	230cm	5/box	000867
Inflation Device					٠
CONMED Balloon Inflation Device				1/box	10085

CONMED Corporation 11311 Concept Blvd Largo, FL 33773 Customer Service: 1-866-4CONMED International Sales: 727-214-3000

Fax: 1-800-438-3051

customerexperience@conmed.com internationalorders@conmed.com