


# SmartNail<sup>®</sup>

FOR REPAIR OF OSTEOCHONDRAL LESIONS  
& FRAGMENT FIXATION

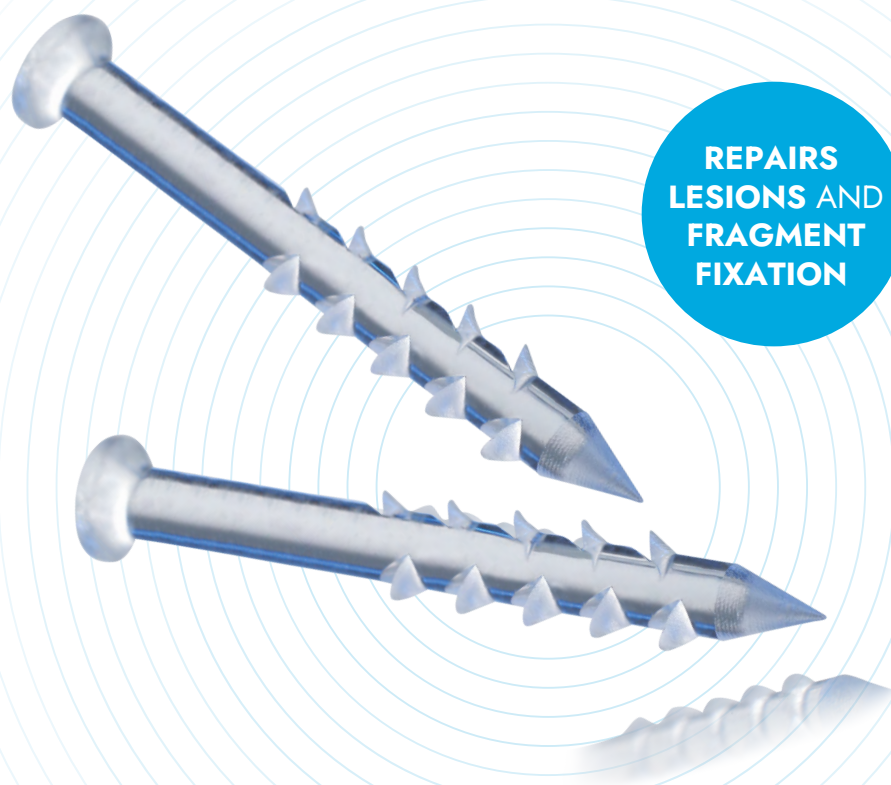


REPAIRS  
LESIONS AND  
FRAGMENT  
FIXATION

**An innovative, clinical solution for  
OCD lesions and fragment fixations  
of non-load bearing bones**

The SmartNail<sup>®</sup> headed implant is designed to offer the compression needed during healing, without requiring a second surgery once healing is complete. The technique offers guidance for simple and quick insertion.

## An Easy Resorbable Solution

- Proprietary Self-Reinforced™ 96L/4D PLA copolymer provides the strongest resorbable implant available
- Effective alternative to metal implants, autogenous transfer systems or regeneration techniques
- Low-profile, round, flat head offers more compression across the fracture when compared to headless implants
- Barbs along the distal portion of the implant are designed to increase compression and pull-out
- 2mm countersink feature integrated in technique
- Available in 1.5mm and 2.4mm diameters

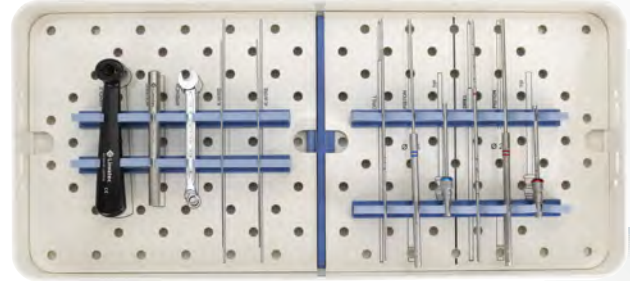

SmartNail® Sterilization Tray



2mm Countersink

## ORDERING INFORMATION

Catalog No.	Product Description	Size
<b>IMPLANTS</b>		
531516	SmartNail® (96L/4D) PLA	1.5mm O.D. x 16mm L
531520	SmartNail® (96L/4D) PLA	1.5mm O.D. x 20mm L
531525	SmartNail® (96L/4D) PLA	1.5mm O.D. x 25mm L
532416	SmartNail® (96L/4D) PLA	2.4mm O.D. x 16mm L
532425	SmartNail® (96L/4D) PLA	2.4mm O.D. x 25mm L
532435	SmartNail® (96L/4D) PLA	2.4mm O.D. x 35mm L
532445	SmartNail® (96L/4D) PLA	2.4mm O.D. x 45mm L

### GENERAL INSTRUMENTS

AA1520-01	Handle
AA1520-04	Housing
AA1520-05	Wrench

Catalog No.	Product Description
<b>1.5MM SET</b>	
SN1520-021	1.5mm Percutaneous Tip (Blue)
SN1520-031	1.5mm Percutaneous Piston (Blue)
SN1520-061	1.5mm Percutaneous Drill Bit (Blue)
<b>2.4MM SET</b>	
PT2400	2.4mm Percutaneous Tip (Red)
PP2400	2.4mm Percutaneous Piston
DB2401	2.4mm Percutaneous Drill Bit

### STERILIZATION TRAY

ST1007	Sterilization Tray (new plates)
--------	---------------------------------