



FX Wire[®]

Advanced Measurement Guidewire

Designed to provide easy access to challenging biliary anatomy, these guidewires combine a clear measurement system with a flexible, 5cm hydrophilic coated radiopaque distal tip for a simple and streamlined clinical experience.



COMPLETECONTROL[™]
Sphincterotome



GORE[®] VIABIL[®]
Short Wire Biliary Endoprosthesis



FIT[®]
Biliary Stone Balloon

To learn more, call
1-866-4CONMED (Toll Free)
727-214-3000 (International)
or visit www.CONMED.com.

FX Wire® Advanced Measurement Guidewire

Benefits Beyond Measurement

When it comes to device delivery in tough-to-reach biliary anatomy, our FX Wire® fully-integrated measurement guidewires provide the size and simplicity necessary for easy instrument access and guidance.

Hydrophilic, Flexible, Atraumatic Tip

Facilitates access past obstructions and tortuous anatomy

REGILIANT™ Nitinol Core


Kink-resistant, optimal column strength provides enhanced stricture access, easy device exchange and superior torque response

Lubricious Silicone Coating

Provides smooth, seamless device exchange

Light Blue PTFE Jacket with Black Markings

Enhances visualization by reducing endoscopic glare

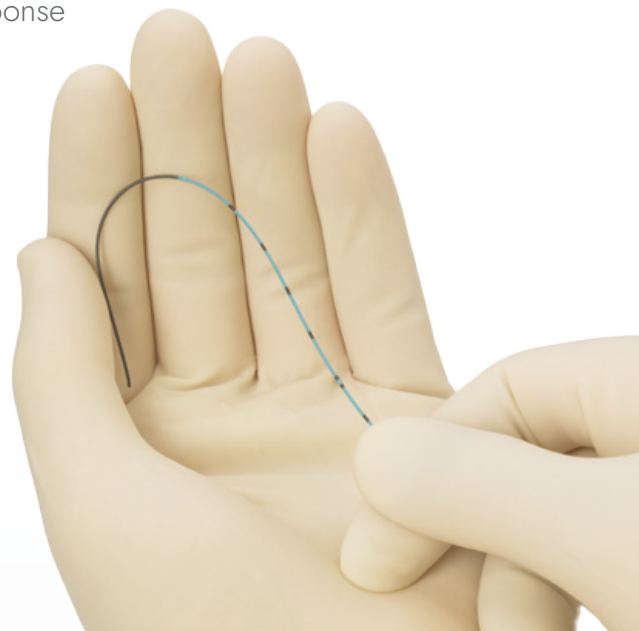

*Hydrophilically coated 5cm
radiopaque floppy tip*

*1cm spaced endoscopic markings
extend 30cm from wire distal tip*

*Proprietary REGILIANT™
Nitinol Core*

PTFE jacket with lubricious coating

*10cm and 15cm Fluoroscopic
marker bands*

*1cm spaced proximal markings
extend from 150cm to 240cm*



Endoscopic Measurement Made Easy

Endoscopic Measurement System

Individual 1cm spaced markings and multiple markings every 5cm for ease of measurement and to help detect inadvertent movement

Radiopaque Marker Bands

Featuring Radiopaque Marker Bands at 10cm and 15cm, which helps facilitate accurate fluoroscopic stricture measurement

To order or learn more, contact your
CONMED Representative or visit
www.CONMED.com.



FX Wire[®]

Advanced Measurement Guidewire

Ordering Information

Description	Guidewire O.D	Length	Tip Type	Torque Vise	Packaging	Catalog Number
Standard Wire	.035" / .89 mm	450 cm	Straight Tip	NO	3/box	FX35450S
Standard Wire	.035" / .89 mm	260 cm	Straight Tip	NO	3/box	FX35260S



CONMED Corporation
11311 Concept Blvd.
Largo, FL 33773

Toll Free: 1-866-4CONMED
International: 727-214-3000
www.CONMED.com

customerexperience@conmed.com
internationalorders@conmed.com