

INFINITY AIM™*

Meniscal Repair Device

Maximize Your Repair with Minimal Anatomic Disruption

When you encounter a meniscal tear, you need tools to efficiently facilitate a repair while preserving anatomic structures. AIM™ provides a strong, low-profile solution using a safe and simple technique.



**AIM™ for a better meniscal repair.
Request a trial for your next case!**

 **CONMED**


* Device is not available for sale in all regions

An all-inside device designed and refined by world renowned knee surgeons. The result? A reliable, easy to use system aimed to perform.


Minimal Anatomic Disruption

In side-by-side testing, AIM™'s hidden implant design reduced the piercing diameter* by up to **36%**¹ when compared to leading competitive devices on the market.


CONMED
AIM™

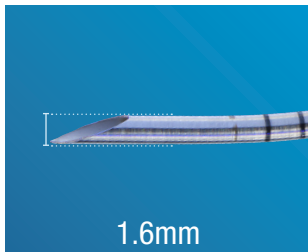


DePuy Synthes
TRUESPAN

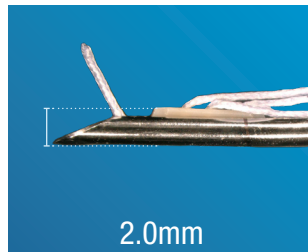

Smith & Nephew
Fast-Fix 360

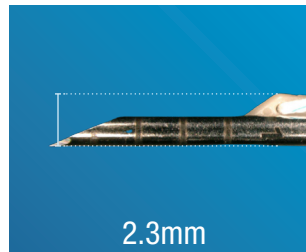

Stryker
AIR+

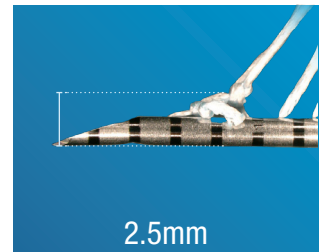


* Piercing diameter defined as the highest diameter that would pierce the meniscus with an exposed implant and needle combined

¹PDD1924378

Designed with Safety in Mind


Retractable Sheath that fully covers the needle allowing for sled-less joint entry without the risk of damaging native anatomy.

Blunt Push Rod extends out past the needle after the second implant is deployed to protect the meniscus while tensioning.

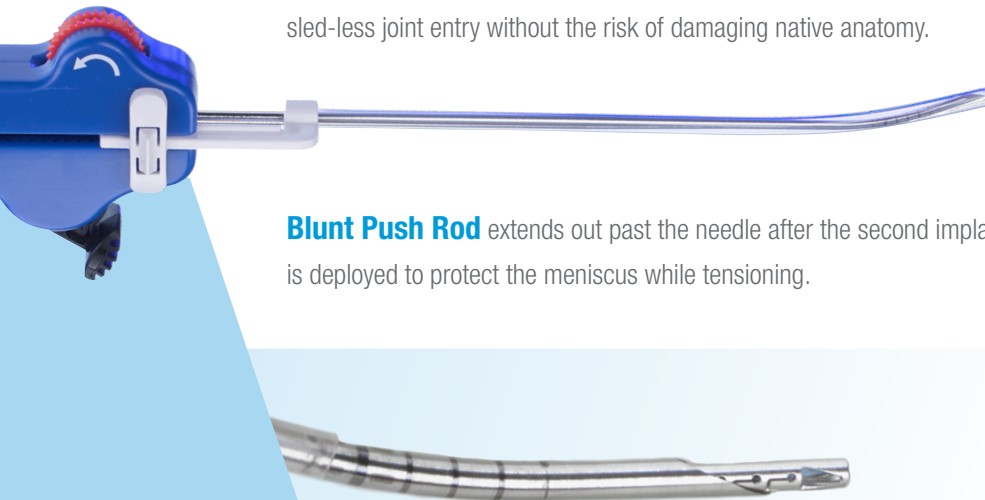

Simplicity

Seamless one-handed deployment technique and in-joint tensioning mechanism promote a simple, reliable repair.

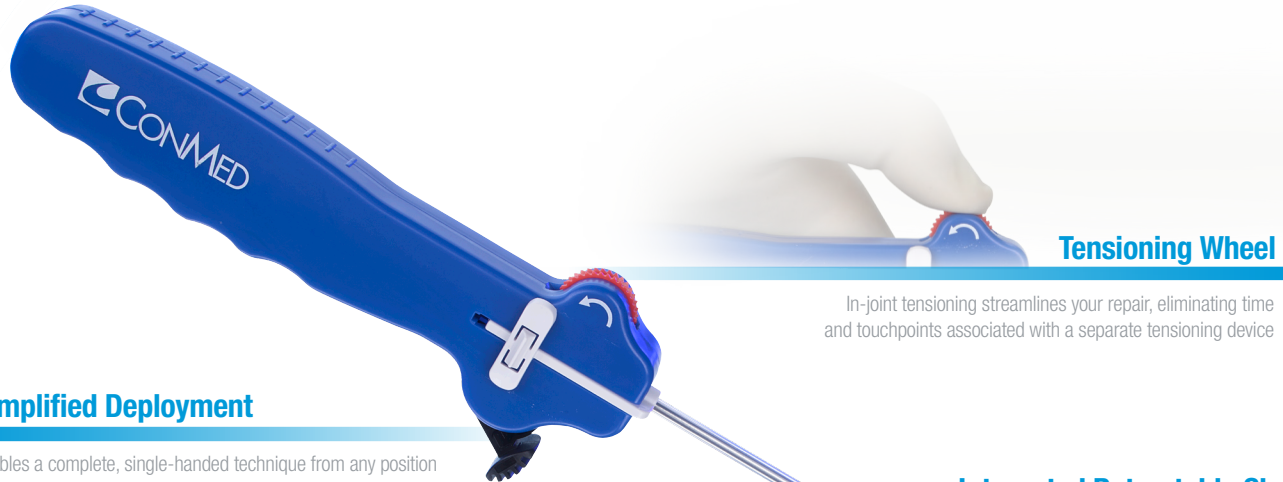

Tensioning Wheel

In-joint tensioning streamlines your repair, eliminating time and touchpoints associated with a separate tensioning device

Simplified Deployment

Enables a complete, single-handed technique from any position




Integrated Retractable Sheath

Easily adjust sheath from full coverage to 20mm of needle exposure for optimal piercing depth

Hidden Implant Within Needle

Preserves the meniscus while piercing

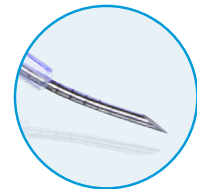

Infinity AIM™ Suture Cutter*



Tapered PEEK Tip

For anatomic protection

* Device is not available for sale in all regions

Infinity AIM™ Meniscal Repair Device

Ordering Information

Description	Catalog Number
Infinity AIM™, 25 Degree Standard Curve	KMR2S
Infinity AIM™, 15 Degree Reverse Curve	KMR2R
Infinity AIM™, Suture Cutter	KSC20

**A Full Offering to Support Your Meniscal Needs.
A to Z, We've Got A Solution For Knee.**



Meniscus Allograft
Transplant system



Zone Specific® II



CONMED Corporation
11311 Concept Blvd.,
Largo, FL 33773

Toll Free: 1-866-4CONMED
International: 727-214-3000
www.CONMED.com

customerexperience@conmed.com
internationalorders@conmed.com