



## LaparoVue® Visibility System Laparoscopic Solutions

Understanding a surgeon's need to see clearly during laparoscopic surgery, the LaparoVue Visibility System is an all-in-one solution designed to take the complexity out of achieving optimal visualization. An easy to use, disposable device that warms, white balances, cleans and defogs laparoscopes.



**AIRSEAL® iFS**  
ACCESS MANAGEMENT SYSTEM



**PLUMEPEN®**  
SMOKE EVACUATION PENCILS



**VISICLEAR®**  
SMOKE EVACUATOR



To learn more about these and other innovative products, call **866-4CONMED** or visit **www.CONMED.com**.

## LaparoVue®

All-in-One Solution for Improved  
Laparoscopic Visibility

### Warms and White Balances


- Warms a laparoscope in 5 minutes to prevent fogging
- Scope Cradle supports and protects during warming and white balancing
- Designed for stability, LaparoVue eliminates surfactant leaks while warming and white balancing

### Cleans and Defogs

- Provides crisp, clear view in 5 seconds
- Offers a chemical and mechanical clean for superior residue removal while helping to prevent intra-operative laparoscope fogging
- Cleaning and defogging port designed to retain surfactant
- Radiopaque VueTip™ Trocar Swabs remove debris and excess liquid from trocar cannula
- Microfiber Cloth captures large residue and excess surfactant



## Distinctive By Design

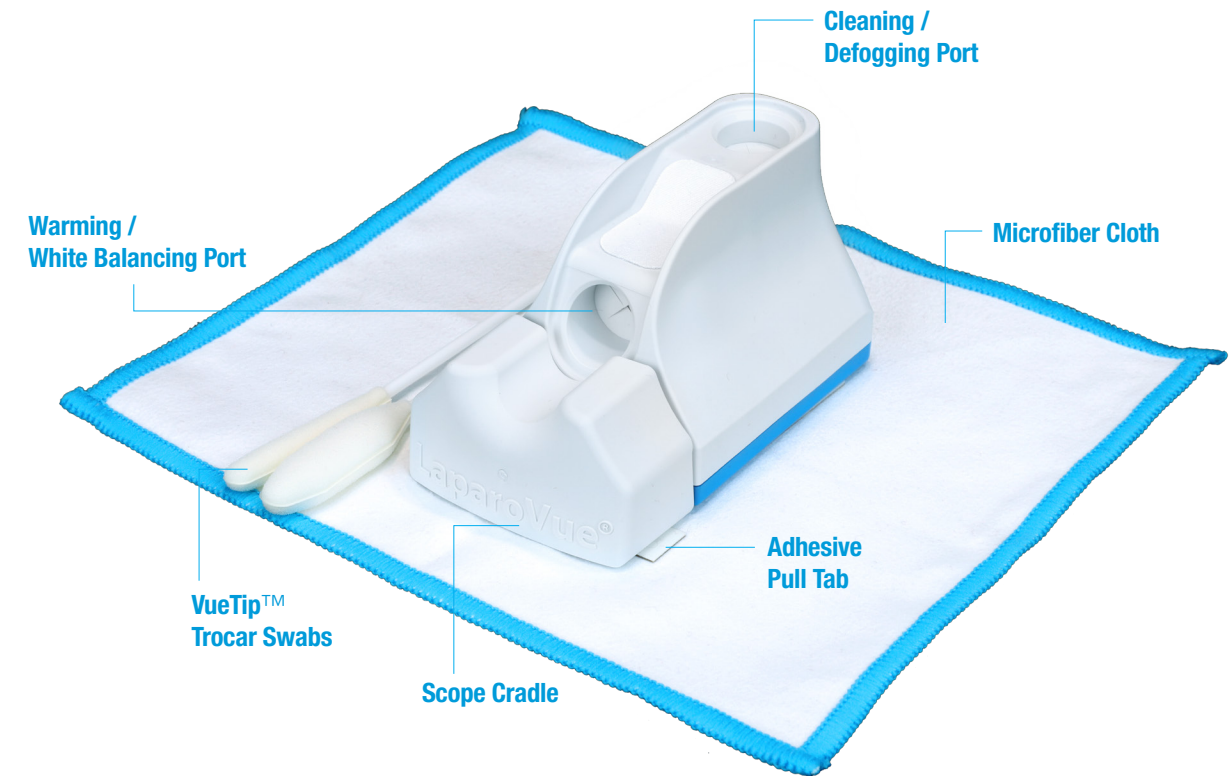

### LaparoVue Product Features

#### One Unit, Two Ports

- Warm a laparoscope while cleaning another, reduces preparation time when changing laparoscopes
- Accepts 3 – 12mm laparoscopes without the need for a port size converter
- Warms, white balances, cleans and defogs for up to 5 hours

#### Designed for Stability

- Scope Cradle safely and securely supports laparoscope to reduce the risk of damage
- Upright functionality reduces chances of surfactant spillage

#### Easy to remove batteries for disposal

#### No assembly requirement streamlines procedure set-up time

#### Single use for maximum infection control



LAPAROVUE®  
VISIBILITY  
SYSTEM

## Ordering Information

Description	Qty	Catalog Number
LaparoVue® Visibility System	10/Box	LAPVUE10
1 LaparoVue® Unit (Main Unit and Scope Cradle)		
2 VueTip™ Trocar Swabs (one large and one small)		
1 Microfiber Cloth		
VueTip™ Trocar Swab Accessory Kit	10/Box	VUETIP10
2 VueTip™ Trocar Swabs (one large and one small)		
1 Microfiber Cloth		

CONMED Corporation  
525 French Road, Utica, NY 13502

866-4CONMED  
International: 727-214-3000

www.CONMED.com  
Customerexperience@CONMED.com