




Paradigm™ ClickLok™* All-Suture Anchor System

Hip Arthroscopists never compromise on performance, which is why we're proud to introduce the ClickLok™ All-Suture Anchor System – CONMED's newest solution for acetabular labral fixation. Developed by some of the world's leading hip surgeons, ClickLok™ is designed to offer the simplicity and reliability you deserve, with a one-handed deployment technique.

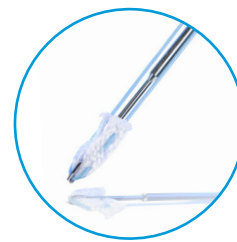

Stay in Control with One-Handed Deployment

The ergonomic delivery system combined with our mechanical deployment technique allows you to deploy the anchor with one hand for a seamless anchor insertion process that eliminates the need to pass control of your scope.

Simple

- Curved delivery system facilitates easy access to the acetabular rim for ideal anchor placement
- Funneled drill guide is designed to help the drill easily and quickly find its insertion point
- Fluted drill bit designed to evacuate bone more efficiently from the tunnel
- Positive stop on the drill bit, when used with the guide, provides simple confirmation that drill depth is achieved

Reliable


The mechanical deployment mechanism ensures consistent and secure fixation, with an audible indicator that confirms final deployment.

Bone Sparing

The 1.8mm ClickLok™ is designed to provide secure fixation with a minimal footprint, helping to preserve crucial bone around the acetabular rim.



Ready to click your way to a better labral repair?
Call 1-866-4CONMED to schedule a demonstration or
visit www.CONMED.com/ParadigmClickLok

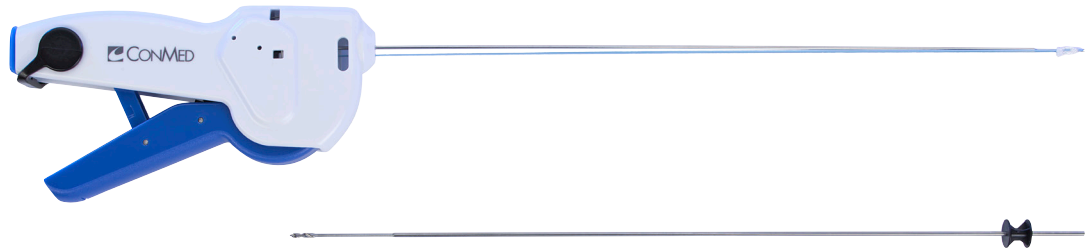

Ordering Information*

Description	Implant / Disposable / Reusable	Catalog Number
Paradigm™ ClickLok™ Anchor, Blue #2 Hi-Fi® Suture	Implant	PASA01
Paradigm™ ClickLok™ Anchor, Black #2 Hi-Fi® Suture	Implant	PASA02
Paradigm™ ClickLok™ Drill Bit	Disposable	PASDB
Paradigm™ ClickLok™ Drill Guide, Straight Crown	Reusable	PASG01A
Paradigm™ ClickLok™ Drill Guide, Straight Fishmouth	Reusable	PASG02A
Paradigm™ ClickLok™ Drill Guide, Curved Crown	Reusable	PASC01A
Paradigm™ ClickLok™ Drill Guide Obturator	Reusable	PASGB01A

*Items not available in all regions

CONMED Corporation
11311 Concept Blvd.,
Largo, FL 33615

Toll Free: 1-866-4CONMED
International: 727-214-3000
www.CONMED.com

customerexperience@conmed.com
internationalorders@conmed.com