



X Wire®

Next Generation Guidewire

CONMED's XWire® Guidewire features proprietary nitinol core composition with a hydrophilically coated floppy distal tip. XWire® provides optimal fluoroscopic visibility and flexibility for effective transition across the elevator of the duodenoscope.

Explore

Tapered Distal 20cm

- Provides optimum flexibility for effective transition across the elevator of the duodenoscope

Hydrophilically Coated Floppy Distal Tip

- Lubricious coating enhances cannulation and stricture access capability


Examine

Radiopaque Floppy Tip

- Provides optimal fluoroscopic visibility
- Optimized flexibility for maneuvering through tortuous angles of the anatomy

Simple Yet Functional Endoscopic Measurement System

- 1cm spaced visual markings, extend 30cm from wire tip
- Multiple markings every 5cm provide a simple, easy to decipher, endoscopic stricture measurement system
- Helps detect inadvertent wire movement
- Ideal color combination for easy visualization

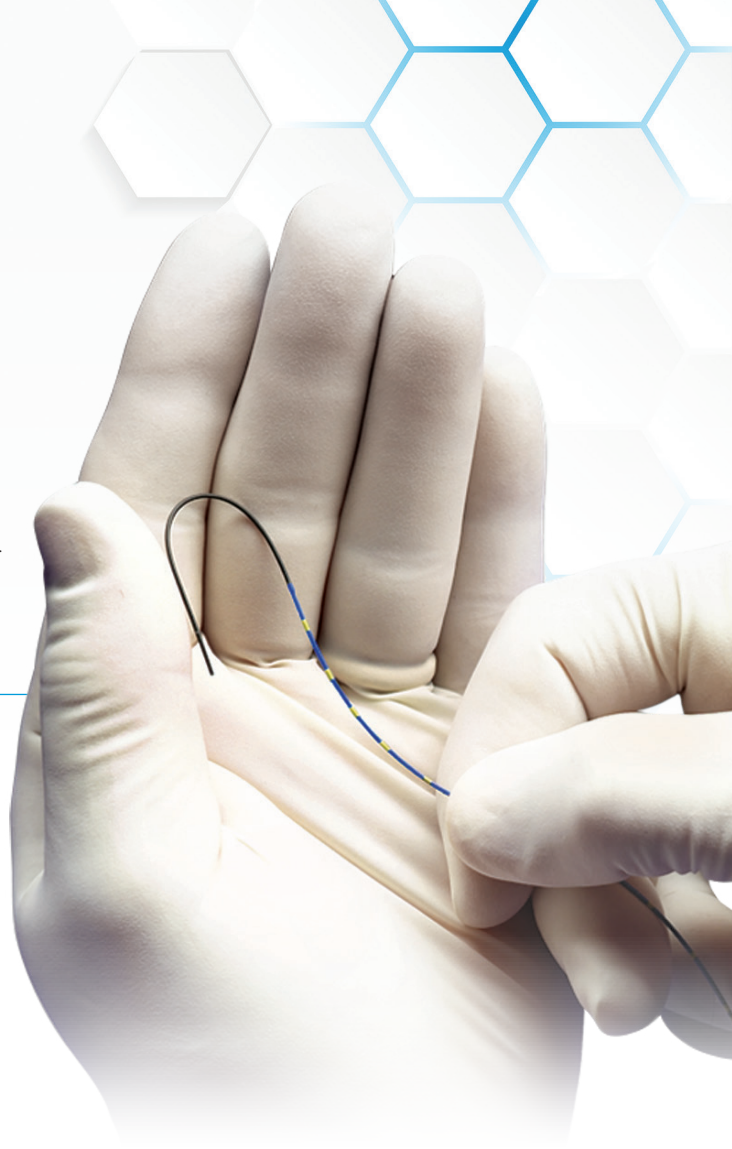

To learn more, call
1-866-4CONMED (Toll Free)
727-214-3000 (International)
or visit www.CONMED.com.





X Wire[®]

Next Generation Guidewire

Ordering Information

Description	O.D	Length	Tip	Compatibility	Vise	Unit	Catalog Number
Standard Wire	.035"	.89mm	450cm	Straight Tip	No	3/box	035450S
Standard Wire	.035"	.89mm	260cm	Straight Tip	No	3/box	035260S
Stiff Wire	.035"	.89mm	450cm	Straight Tip	No	3/box	035450SS
Stiff Wire	.035"	.89mm	260cm	Straight Tip	No	3/box	035260SS
Standard Wire	.025"	.64mm	450cm	Straight Tip	No	3/box	025450S
Standard Wire	.025"	.64mm	260cm	Straight Tip	No	3/box	025260S

CONMED Corporation
11311 Concept Blvd.
Largo, FL 33773

Toll Free: 1-866-4CONMED
International: 727-214-3000
www.CONMED.com

customerexperience@conmed.com
internationalorders@conmed.com