



Endoscope Cleaning Brushes

Engineered for optimal cleaning and convenience, ConMed has a broad portfolio of disposable and reusable endoscope cleaning brushes that feature a universal design to meet your clinical needs.



SCOPE SAVER®
BITE BLOCK



CLEARGUARD™
BIOPSY VALVE



OPTIMIZER®
MULTI-CHAMBER POLYP TRAP



To learn more about these and other innovative products, call 800-448-6506 or visit ConMed.com.

ENDOSCOPE CLEANING BRUSHES




ROUNDED BRUSH TIP

Scrub Smarter.

ConMed's Endoscope Cleaning Brushes feature a universal design that's been created to effectively and thoroughly clean a full range of endoscope channels – helping to ensure proper cleaning. All of the brushes in this comprehensive portfolio feature a rounded brush tip to facilitate passage through the scope channel while protecting the biopsy channel from damage.

Channel Master™ Disposable Endoscope Cleaning Brushes

Our Channel Master™ series of disposable endoscope cleaning brushes feature a custom crimped shaft designed to effectively traverse tight channel bends and deflect through 90° valve port transitions. To help ensure proper cleaning of the valves, button and biopsy channels, ConMed offers a combination brush featuring a larger head of soft bristles with a convenient “No Slip” handle at one end and channel brush at the other. While the dual ended brush utilizes two brushes with one catheter pass.

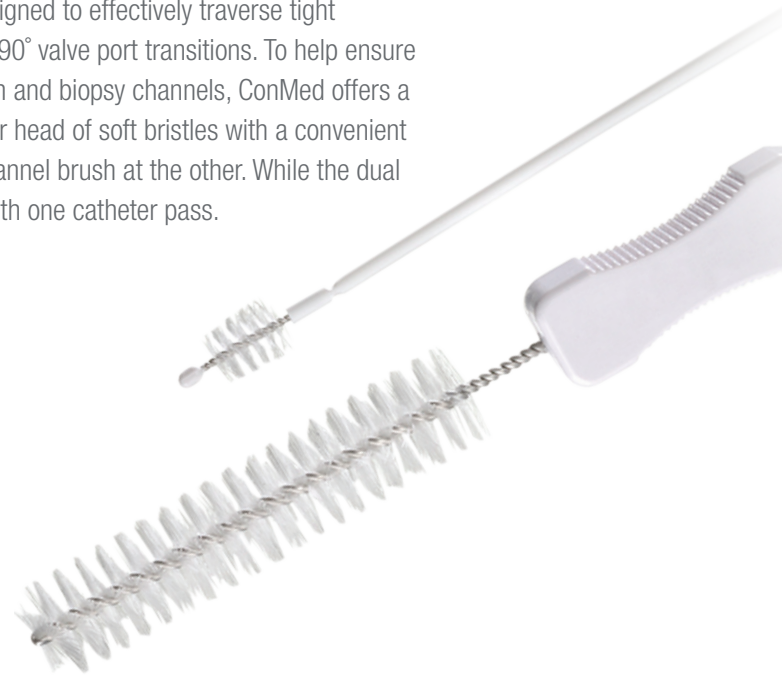



Blue Bullet® Disposable Cleaning Brushes

Our complete Blue Bullet® series of disposable cleaning brushes feature blue opaque tubing, a longer working shaft and stiffer bristles to effectively clean all endoscope channels. ConMed offers a combination Blue Bullet Brush which features a 38mm head with a thumb hold handle to help ensure proper cleaning of the larger valves/channels and 13mm channel brush at the other end.

Whistle® Disposable Endoscope Cleaning Brush

Our Whistle® disposable endoscope cleaning brush features a female luer adaptor, hollow shaft and port hole for fluid injection to effectively clean channel obstacles.

Reusable Endoscope Cleaning Brushes

Our Reusable Cleaning Brushes feature a controllable ring handle and a Teflon® coated spring sheath designed for multiple uses.

ENDOSCOPE
CLEANING
BRUSHES



WHISTLE® BRUSH



REUSABLE BRUSHES

Ordering Information

Scope	Brush Diameter Length		Shaft Length	Unit	Catalog Number
Channel Master™ Endoscope Cleaning Brushes					
Non-Sterile, Disposable – Biopsy Channel Range 2.0-4.2mm, Shaft O.D. 1.8mm					
Pediatric Brush	6.4mm	7.4mm	160cm	50/box	000610
Dual-Ended Brush	6.4mm	7.4mm	230cm	50/box	000615
Single-Ended Brush	6.4mm	7.4mm	230cm	50/box	000639
Combination Valve/Cleaning Brush for Gastroscope/Colonoscope	7.9mm/ 6.4mm	29.2mm/ 7.4mm	230cm	50/box	000617
Large Combination Valve/Cleaning Brush for Gastroscope/Colonoscope	9.7mm/ 6.4mm	45.2mm/ 7.4 mm	230cm	50/box	000619*
* Large combination brush includes longer & larger diameter bristles on valve brush.					
Whistle® Cleaning Brush					
Non-Sterile, Disposable – Biopsy Channel Minimum 2.8mm, Shaft O.D. 1.8mm					
Whistle® Cleaning Brush	6.4mm	13.1mm	230cm	50/box	000101
Contains female luer adaptor, hollow shaft, and port hole for fluid injection to soften channel obstructions.					
Blue Bullet® Disposable Cleaning Brushes					
Non-Sterile, Disposable – Biopsy Channel Range 2.0-5.5mm					
Bronchoscope	3mm	13mm	115cm	50/box	480
Gastro/Duodeno/Colonoscope	3mm	13mm	250cm	50/box	473
Gastro/Duodeno/Colonoscope	4mm	13mm	250cm	50/box	474
Gastro/Duodeno/Colonoscope	5mm	13mm	250cm	50/box	478
Gastro/Duodeno/Colonoscope	6mm	13mm	250cm	50/box	476
Double End Cleaning Brush for Gastroscope/Colonoscope	5mm	13mm	250cm	50/box	478D
Combination Valve/Cleaning Brush for Gastroscope/Colonoscope	10mm/ 5mm	38mm/ 13mm	250cm	50/box	478C
The 400 series Blue Bullet® feature a blue opaque tubing.					
Reusable Cleaning Brush					
Non-Sterile, Disposable					
Bronchoscope	3mm	9.5mm	115cm	1/box	380
Gastro/Duodeno/Colonoscope	6mm	9.5mm	230cm	1/box	378
The reusable scope cleaning brush features a controllable ring handle and a Teflon® coated spring sheath..					

Customer Service: 1-800-448-6506
International Sales: 1-315-797-8375
Fax: 1-800-438-3051

ConMed.com
info@mail.conmed.com